

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

# **Product Sheet**

## **H\_CDH17**(ΔEC1,Flag-EC2-7) HEK-293 Cell Line

Catalog number: GM-C31203

Version 3.3.1.250808

 $\label{eq:Line} \mbox{$H_{CDH17(\Delta EC1,Flag-EC2-7)$ HEK-293$ Cell Line is a clonal stable HEK-293$ cell line} \mbox{$Description}$ 

constitutively expressing human CDH17.

**Quantity** 5E6 Cells per vial,1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Human CDH17(ΔEC1,Flag-EC2-7) & N-Flag

**Gene ID/Uniprot ID** Q12864( AA Ser 129- Gly 777; AA Ala 788 - Ile 808)

Host Cell HEK-293

**Recovery Medium** DMEM +10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+0.75 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

**Growth properties** Adherent

Growth Conditions 37°C, 5% CO<sub>2</sub>

**Mycoplasma Testing** The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Materials**

Reagent	Manufacturer/Catalogue No.
Puromycin	Genomeditech/GM-040401
Pen/Strep	Thermo/15140-122
Fetal Bovine Serum	ExCell/FSP500
DMEM	Gibco/C11995500BT
Anti-Flag mIgG1 Antibody	Genomeditech/GM-30726AB

### **Figures**

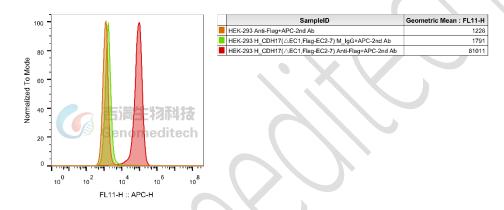


Figure 1 |  $H_CDH17(\Delta EC1, Flag-EC2-7)$  HEK-293 Cell Line (Cat. GM-C31203) was determined by flow cytometry using Anti-Flag mIgG1 Antibody (Cat. GM-30726AB).

#### **Cell Recovery**

Recovery Medium: DMEM +10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

## **Cell Freezing**

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

## Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

#### **Notes**

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Sequence**

Flag-CDH17 Q12864( AA Ser 129- Gly 777; AA Ala 788 - Ile 808)

DYKDDDDKGSSKYEGSVRQNSRPGKPFLYVNATDLDDPATPNGQLYYQIVIQLPMINNVMYFQINNKTGAIS
LTREGSQELNPAKNPSYNLVISVKDMGGQSENSFSDTTSVDIIVTENIWKAPKPVEMVENSTDPHPIKITQVRW
NDPGAQYSLVDKEKLPRFPFSIDQEGDIYVTQPLDREEKDAYVFYAVAKDEYGKPLSYPLEIHVKVKDINDNP
PTCPSPVTVFEVQENERLGNSIGTLTAHDRDEENTANSFLNYRIVEQTPKLPMDGLFLIQTYAGMLQLAKQSL
KKQDTPQYNLTIEVSDKDFKTLCFVQINVIDINDQIPIFEKSDYGNLTLAEDTNIGSTILTIQATDADEPFTGSSK
ILYHIIKGDSEGRLGVDTDPHTNTGYVIIKKPLDFETAAVSNIVFKAENPEPLVFGVKYNASSFAKFTLIVTDVN
EAPQFSQHVFQAKVSEDVAIGTKVGNVTAKDPEGLDISYSLRGDTRGWLKIDHVTGEIFSVAPLDREAGSPYR
VQVVATEVGGSSLSSVSEFHLILMDVNDNPPRLAKDYTGLFFCHPLSAPGSLIFEATDDDQHLFRGPHFTFSLG
SGSLQNDWEVSKINGTHARLSTRHTEFEEREYVVLIRINDGGRPPLEGIVSLPVTFCSCVEGSCFRPAG
AVGILLTTLLVIGIILAVVFI

#### **Related Products**

CDH3		
Cynomolgus_CDH3 CHO-K1 Cell Line	H_CDH3 CHO-K1 Cell Line	
H_CDH3 HEK-293 Cell Line		
Anti-H_CDH3 hIgG1 Antibody		
CDH6		
Cynomolgus_CDH6 CHO-K1 Cell Line	H_CDH6 CHO-K1 Cell Line	
H_CDH6 HEK-293 Cell Line		
Anti-CDH6 hIgG1 Reference Antibody (Ralubio)	Anti-H_CDH6 hIgG1 Antibody(H01L02)	
CDH17		
Cynomolgus_CDH17 HEK-293 Cell Line	Cynomolgus_CDH17(XP_005563762.1) HEK-293 Cell Line	
H_CDH17 CHO-K1 Cell Line	H_CDH17 CT26 Cell Line	
H_CDH17 HCT116 Cell Line	H_CDH17 HEK-293 Cell Line	
H_CDH17 LLC1 Cell Line	H_CDH17 MC38 Cell Line	
H_CDH17 RKO Cell Line	H_CDH17 SW480 Cell Line	
H_CDH17(ΔEC1-2,Flag-EC3-7) HEK-293 Cell Line	H_CDH17(ΔEC1-3,Flag-EC4-7) HEK-293 Cell Line	
H_CDH17(ΔEC1-4,Flag-EC5-7) HEK-293 Cell Line	H_CDH17(ΔEC1-5,Flag-EC6-7) HEK-293 Cell Line	
H_CDH17(ΔEC1-6,Flag-EC7) HEK-293 Cell Line	Mouse_CDH17 HEK-293 Cell Line	
Rat_CDH17 HEK-293 Cell Line	Rhesus_CDH17 HEK-293 Cell Line	
Anti-CDH17 hIgG1 Antibody(BI-905711)	Anti-CDH17 hIgG1 Antibody(VHHI-28BB)	
Anti-CDH17 hIgG1 Reference Antibody(BI-905711)		
Biotinylated Human CDH17 Protein; His-Avi Tag	Cynomolgus CDH17 Protein; His Tag	
Human CDH17 Protein; His Tag	Mouse CDH17 Protein; His Tag	



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

#### **License Agreement:**

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope;
   modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.