

Product Sheet

Membrane Bound H_IL18 CHO-K1 Cell Line

Catalog number: GM-C44216

Version 3.3.1.260602

Description	Membrane Bound H_IL18 CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively expresses the human IL18 gene, constructed using lentiviral technology.
Quantity	5E6 Cells per vial, 1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt
Target	Human_IL18
Gene ID/Uniprot ID	Q14116-1(AA Tyr 37 - Asp 193)
Host Cell	CHO-K1
Recovery Medium	F12K+10% FBS+1% P.S
Growth medium	F12K+10% FBS+1% P.S+4 µg/mL Puromycin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-IL18 hIgG1 Antibody(CMK389)	Genomeditech/GM-88326AB
H_IL18 Reporter 293 Cell Line	Genomeditech/GM-C21147

Figures

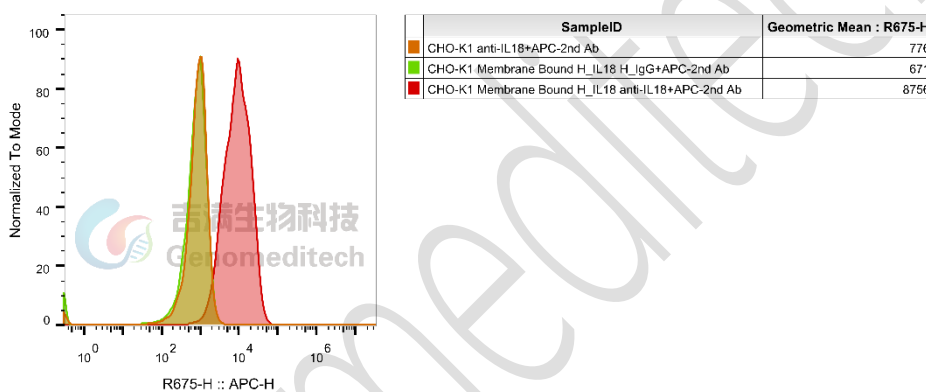


Figure 1 | Membrane Bound H_IL18 CHO-K1 Cell Line (Cat. GM-C44216) was determined by flow cytometry using Anti-IL18 hIgG1 Antibody(CMK389) (Cat. GM-88326AB).

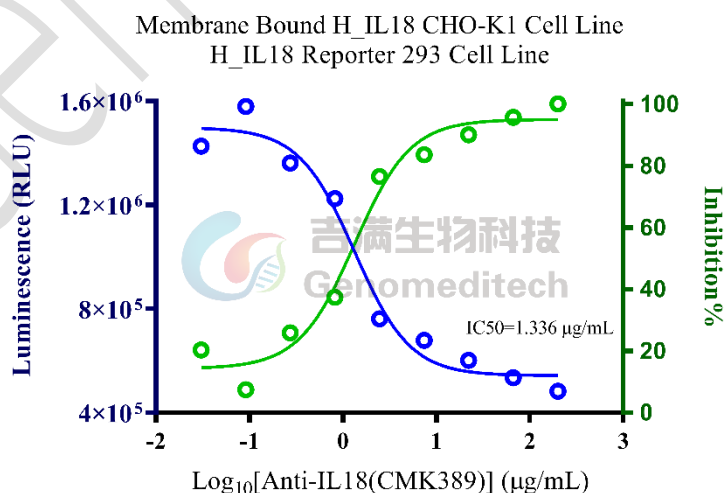


Figure 2 | Response to Anti-IL18 hIgG1 Antibody(CMK389). H_IL18 Reporter 293 Cell Line (Cat. GM-C21147) was seeded at a density of 1.5E4 cells per well in a 96-well plate and incubated overnight. The next day, serial dilutions

of the Anti-IL18 hIgG1 Antibody(CMK389) (Cat. GM-88326AB) were incubated with 1E3 cells/well of the Membrane Bound H_IL18 CHO-K1 Cell Line (Cat. GM-C44216) in a 96-well plate for 10 min, and then added to the pre-seeded cells. The mixture was incubated for an additional 16 hours. Firefly luciferase activity was measured using the Luciferase Reporter Assay Kit (Genomeditech) (left Y-axis, relative luminescence units), with inhibition percentages shown on the right Y-axis. Data are shown by drug mass concentration.

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at 176 x g for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.

- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

IL18 Q14116-1(p.37-193)

YFGKLESKLSVIRNLNDQVLFIDQGNRPLFEDMTDSDCRDNAPRTIFIISMYKDSQPRGMAVTISVKCEKISTLS
CENKIISFKEMNPPDNIKDTKSDIIFQRSVPGHDNKMQFESSYEGYFLACEKERDLFKLILKKEDELGDRSIM
FTVQNE

Related Products

IL-18	
H_IL18 Reporter 293 Cell Line	Mouse_IL18 Reporter 293 Cell Line
Anti-IL18BP hIgG1 Antibody(SA338a)	Anti-IL18RA mIgG1 Antibody(1G12)
Anti-IL18RB hIgG1 Antibody (P20-V1)	
Biotinylated Human IL-18BP Protein; His-Avi Tag	Biotinylated Human IL-18R beta Protein; hFc-Avi Tag
Human IL-18BP Protein; hFc Tag	Human IL-18R alpha Protein; hFc Tag
Human IL-18R beta Protein; hFc Tag	

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly used in experiments involving humans.

- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope; modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.

Genomeditech