

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

H_GDF15 Reporter 293 Cell Line


Catalog number: GM-C06718

Version 3.3.1.250124

Growth differentiation factor 15 (GDF15) is a distant member of the transforming growth factor (TGF)- β family and is a secreted protein that circulates as a 25 kDa dimer. It is a protein belonging to the transforming growth factor β superfamily.

The glial cell line-derived neurotrophic factor (GDNF) family receptor alpha (GFRAL) is the sole receptor for GDF15. GDF15 binds to GFRAL with high affinity but does not interact with other members of the GDNF receptor family. Upon binding GDF15, GFRAL forms a complex with the co-receptor RET, enabling the observation of typical GDNF family cell signaling, which is associated with the phosphorylation of RET, AKT, ERK, and PLC- γ 1.

H_GDF15 Reporter 293 Cell Line is a clonal stable HEK-293 cell line constructed using lentiviral technology, constitutive expression of the human GFRAL gene and human RET, along with signal-dependent expression of a luciferase reporter gene. When GDF15 binds to GFRAL-RET, it activates downstream signaling pathways, leading to the expression of luciferase. The luciferase activity measurement indicates the activation level of the signaling pathway and can thus be used to evaluate the in vitro effects of drugs related to this signaling pathway.

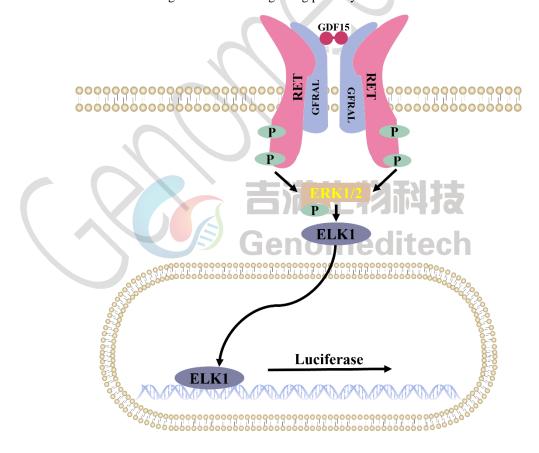



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Specifications

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Recovery Medium DMEM+10% FBS+1% P.S

Hygromycin+0.75 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/GM-040404
G418	Genomeditech/GM-040402
Hygromycin	Genomeditech/GM-040403
Puromycin	Genomeditech/GM-040401
Human GDF15 Protein; His Tag	Genomeditech/GM-87627RP
GMOne-Step Luciferase Reporter Gene Assay Kit	Genomeditech/GM-040503

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Figures

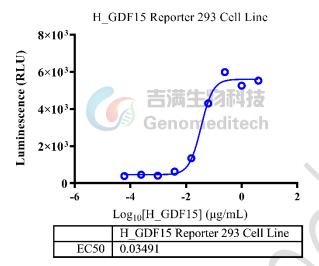


Figure 1 | Response to Human GDF15 Protein; His Tag. The H_GDF15 Reporter 293 Cell Line (Cat. GM-C06718) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Human GDF15 Protein; His Tag (Cat. GM-87627RP) in assay buffer (DMEM+1% FBS+1% P.S) for 16 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). The maximum induction fold was approximately [14.0]. Data are shown by drug mass concentration.

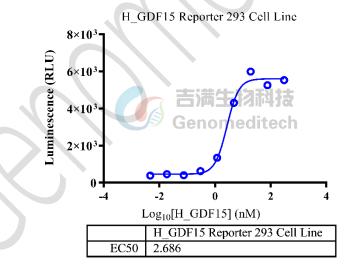


Figure 2 | Response to Human GDF15 Protein; His Tag. The H_GDF15 Reporter 293 Cell Line (Cat. GM-C06718) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Human GDF15 Protein; His Tag (Cat. GM-87627RP) in assay buffer (DMEM+1% FBS+1% P.S) for 16 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). The maximum induction fold was approximately [14.0]. Data are shown by drug molar concentration.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C.

Storage at -70°C will result in loss of viability.

Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring

and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).

b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying

with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.

c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately

176 x g for 5 minutes. Discard supernatant.

d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium

described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

Centrifuge at 176 x g for 3 minutes to collect cells.

Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL. b)

Aliquot 1 mL into each vial. c)

Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid d)

nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+1.5 μg/mL Blasticidin+400 μg/mL G418+125 μg/mL Hygromycin+0.75

μg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth

medium.

Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of

1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability

due to compression.

b) Remove and discard culture medium.

c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.

Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell d)

layer is dispersed (usually within 30 to 60 seconds at 37°C).

Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.

Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.

b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

Related Products

GDF15 GFRAL	
H_GDF15 Reporter 293 DDX35TM Cell Line	Mouse_GDF15 Reporter 293 Cell Line
Membrane Bound H_GDF15 CHO-K1 Cell Line	
Anti-GDF15 hIgG1 Reference Antibody(Ponsbio)	Anti-H_GDF15 hIgG4 Antibody(CTL-002)
Anti-H_GFRAL hIgG4 Antibody(QUEL-0301)	
Biotinylated Human GFRAL Protein; His-Avi Tag	Human GDF15 Protein; His Tag
Human GFRAL Protein; His Tag	Human pro-GDF15 Protein; His Tag

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.