

Product Sheet

H_NKG2A H_CD94 CHO-K1 Cell Line

Catalog number: GM-C15621

Version 3.3.1.241206

Description	H_NKG2A H_CD94 CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively expresses the human NKG2A and human CD94 genes, constructed using lentiviral technology.
Quantity	5E6 Cells per vial, 1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt
Target	Human_NKG2A & Human_CD94
Gene ID/Uniprot ID	NP_998823.2 & NP_001107868.2
Host Cell	CHO-K1
Recovery Medium	F12K+10% FBS+1% P.S
Growth medium	F12K+10% FBS+1% P.S+4 µg/mL Blasticidin+4 µg/mL Puromycin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/ GM-040404
Puromycin	Genomeditech/ GM-040401
Purified anti-human CD94 Antibody	BioLegend/305502
Anti-H_KLRC1 hIgG4 Antibody(Monalizumab)	Genomeditech/ GM-34467AB

Figures

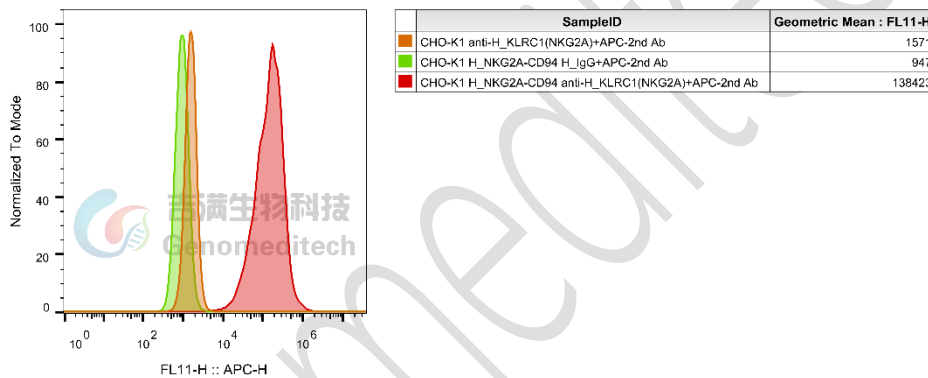


Figure 1 | H_NKG2A H_CD94 CHO-K1 Cell Line (Cat. GM-C15621) was determined by flow cytometry using Anti-H_KLRC1 hIgG4 Antibody (Cat. [GM-34467AB](#)).

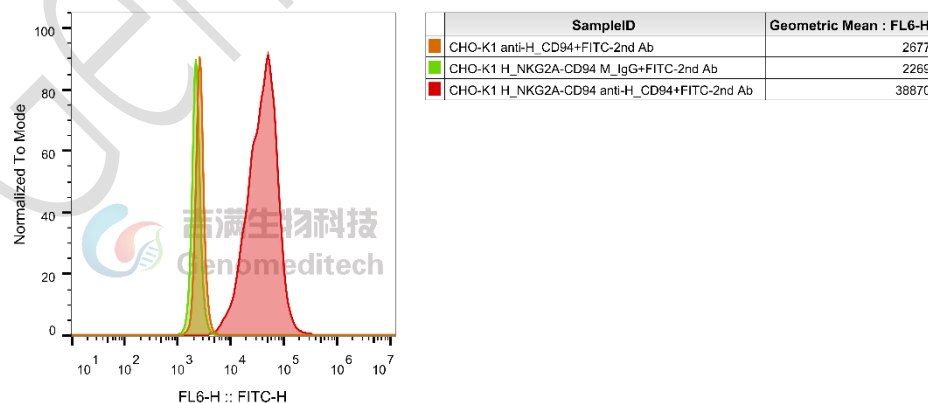


Figure 2 | H_NKG2A H_CD94 CHO-K1 Cell Line (Cat. GM-C15621) was determined by flow cytometry using Purified anti-human CD94 Antibody (BioLegend/305502).

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C . Storage at -70°C will result in loss of viability.

- Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately $176 \times g$ for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- Incubate the culture at 37°C in a suitable incubator. A 5% CO_2 in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at $176 \times g$ for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5×10^6 cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 $\mu\text{g}/\text{mL}$ Blasticidin+4 $\mu\text{g}/\text{mL}$ Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- Incubate cultures at 37°C .

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

KLRC1(NKG2A) [NP_998823.2](#)

MDNQGVIYSDLNLPNPKRQQRKPKGNKNSILATEQEITYAELNLQKASQDFQGNKTYHCKDLPSAPEKLI
 VGILGIICLILMASVVTIVVIPSTLIQRHNSSLNTRTQKARHCGHCPEEWITYNSCYIYGKERRTWEESLLAC
 TSKNSSLLSIDNEEEMKFLSIISPSSWIGVFRNSSHHPWVTMNGLAFLKHEIKSDNAELNCAVLQVNRKLSAQ
 CGSSIIYHCKHKL*

KLRD1(CD94) [NP_001107868.2](#)

MAVFKTTLWRLISGTLGIICLSLMTLGIILLKNSFTKLSIEPAFTPGPNIELQKSDCCSCQEKWVGYRCNCYFI
 SSEQKTWNESRHLCASQKSSLLQLQNTDELDFMSSSQFYWIGLSYSEEHTAWLWENGALSQYLFPSFETF
 NTKNCIAYNPNGNALDESCEDKNRYICKQQLI*

Related Products

NKG2A:HLA-E	
H_HLA-E HEK-293 Cell Line	
Anti-H_KLRC1 hIgG4 Antibody(Monalizumab)	Anti-HLA-E hIgG1 Antibody(ABX-0020)
NKG2D	
H_NKG2D Blockade Reporter Jurkat Cell Line	Cynomolgus_NKG2D CHO-K1 Cell Line
H_NKG2D CHO-K1 Cell Line	H_NKG2D HEK-293 Cell Line
Anti-H_KLRK1(NKG2D) hIgG4 Antibody(Tesnatilimab)	Anti-NKG2D hIgG1 Antibody(A49MI)

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.