

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

# **Product Sheet**

# H HLA-G1 CHO-K1 Cell Line

Catalog number: GM-C15840

Version 3.3.1.241206

H\_HLA-G1 CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively

**Description** expresses the human HLA-G1 and human B2M genes, constructed using lentiviral

technology.

**Quantity** 5E6 Cells per vial,1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Human\_HLA-G1 & Human\_B2M

**Gene ID/Uniprot ID** P17693-1 & P61769

Host Cell CHO-K1

**Recovery Medium** F12K+10% FBS+1% P.S

**Growth medium** F12K+10% FBS+1% P.S+100 μg/mL Hygromycin+4 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

**Growth properties** Adherent

**Growth Conditions** 37°C, 5% CO<sub>2</sub>

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

Note

It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Materials**

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Hygromycin	Genomeditech/GM-040403
Puromycin	Genomeditech/GM-040401
Anti-H_HLA-G1 hIgG1 Antibody(38373)	Genomeditech/GM-28208AB

## **Figures**

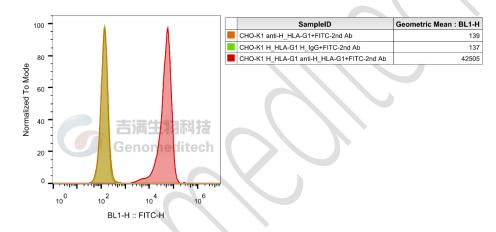


Figure 1 | H\_HLA-G1 CHO-K1 Cell Line (Cat. GM-C15840) was determined by flow cytometry using Anti-H\_HLA-G1 hIgG1 Antibody(38373) (Cat. GM-28208AB).

## **Cell Recovery**

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.

2 / 4



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

# **Cell Freezing**

Freezing Medium: 90% FBS+10% DMSO

a) Centrifuge at 176 x g for 3 minutes to collect cells.

b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.

c) Aliquot 1 mL into each vial.

d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

## Cell passage

Growth medium: F12K+10% FBS+1% P.S+100 µg/mL Hygromycin+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

a) Remove and discard culture medium.

b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.

c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).

d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.

e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.

f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.

g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

#### Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

## **Sequence**

HLA-G1(HLA-G) P17693-1

MVVMAPRTLFLLLSGALTLTETWAGSHSMRYFSAAVSRPGRGEPRFIAMGYVDDTQFVRFDSDSACPRMEP RAPWVEQEGPEYWEEETRNTKAHAQTDRMNLQTLRGYYNQSEASSHTLQWMIGCDLGSDGRLLRGYEQY AYDGKDYLALNEDLRSWTAADTAAQISKRKCEAANVAEQRRAYLEGTCVEWLHRYLENGKEMLQRADPP



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

KTHVTHHPVFDYEATLRCWALGFYPAEIILTWQRDGEDQTQDVELVETRPAGDGTFQKWAAVVVPSGEEQR YTCHVQHEGLPEPLMLRWKQSSLPTIPIMGIVAGLVVLAAVVTGAAVAAVLWRKKSSD

B2M P61769

MSRSVALAVLALLSLSGLEAIQRTPKIQVYSRHPAENGKSNFLNCYVSGFHPSDIEVDLLKNGERIEKVEHSDL SFSKDWSFYLLYYTEFTPTEKDEYACRVNHVTLSQPKIVKWDRDM

#### **Related Products**

HLA-G	
H_HLA-G1 HEK-293 Cell Line	H_HLA-G1 MC38 Cell Line
H_HLA-G1 OKT3 CHO-K1 Cell Line	H_HLA-G1 PDL1 MC38(mouse_PDL1 KO) Cell Line
H_HLA-G1 SK-OV3 Cell Line	Rhesus_MAMU-AG HEK-293 Cell Line
Anti-H_HLA-G1 hIgG1 Antibody(38373)	Anti-H_HLA-G1 mIgG2a Antibody

# **Limited Use License Agreement**

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.