

Product Sheet

H_FCGR2B(CD32B) CHO-K1 Cell Line

Catalog number: GM-C16925

Version 3.3.1.250123

Description	H_FCGR2B(CD32B) CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively expresses the human FCGR2B(CD32B) gene, constructed using lentiviral technology.
Quantity	5E6 Cells per vial, 1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt
Target	Human_FCGR2B(CD32B)
Gene ID/Uniprot ID	P31994-1
Host Cell	CHO-K1
Recovery Medium	F12K+10% FBS+1% P.S
Growth medium	F12K+10% FBS+1% P.S+4 µg/mL Puromycin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
APC anti-human CD32B/C	Biolegend/398304
Anti-H_BDCA2 hIgG1 Antibody(Litifilimab)	Genomeditech/GM-31294AB
H_BDCA2 Reporter Jurkat Cell Line	Genomeditech/GM-C13225
GMOne-Step Luciferase Reporter Gene Assay Kit	Genomeditech/GM-040503

Figures

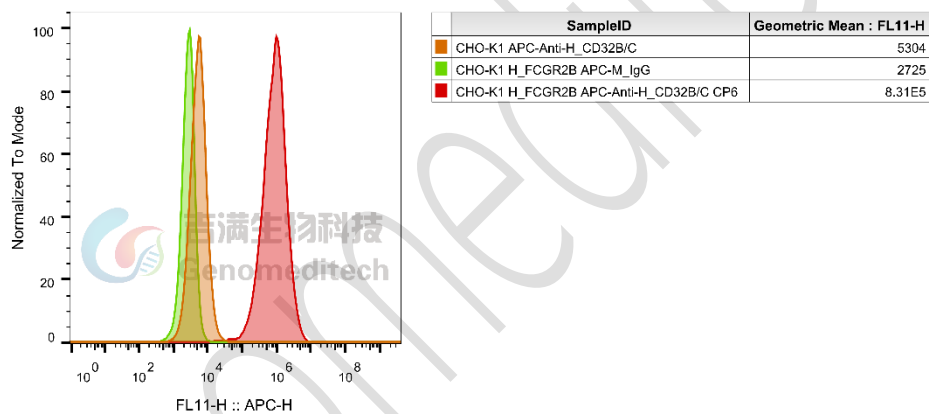


Figure 1 | H_FCGR2B(CD32B) CHO-K1 Cell Line (Cat. GM-C16925) was determined by flow cytometry using APC anti-human CD32B/C (BioLegend/398304).

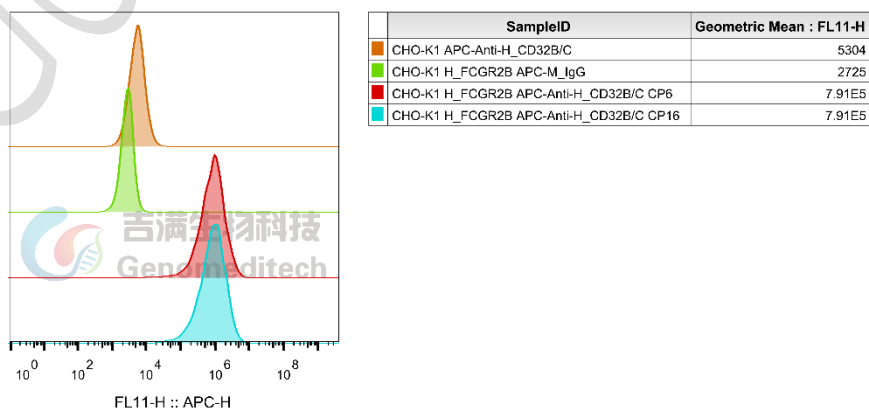


Figure 2 | The passage stability of the H_FCGR2B(CD32B) CHO-K1 Cell Line (Cat. GM-C16925) was determined by flow cytometry using APC anti-human CD32B/C (BioLegend/398304).

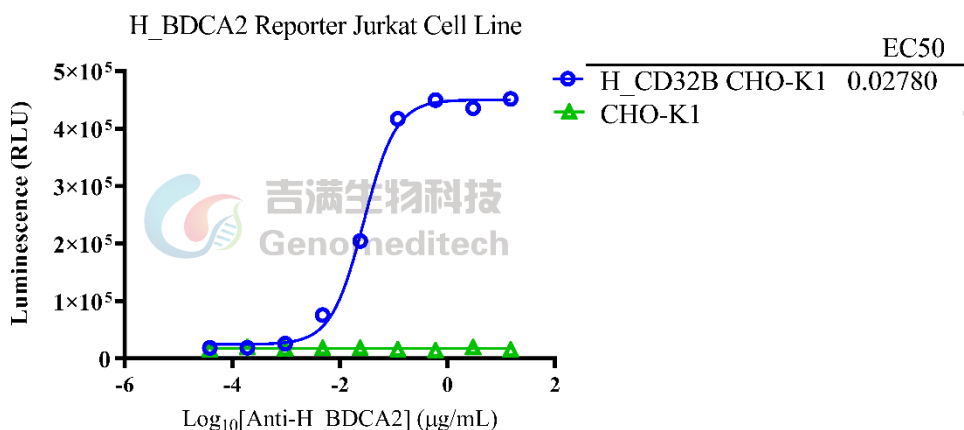


Figure 3 | Response to Anti-H_BDCA2 hIgG1 Antibody. The H_FCGR2B (CD32B) CHO-K1 Cell Line (Cat. GM-C16925) and the CHO-K1 Cell Line were seeded at a density of 1E4 cells per well in a 96-well plate and incubated overnight. The next day, serial dilutions of Anti-H_BDCA2 hIgG1 Antibody (Litifilimab, Cat. GM-31294AB) and the H_BDCA2 Reporter Jurkat Cell Line (Cat. GM-C13225) at a concentration of 1E5 cells per well were added to the pre-seeded cells. The mixture was incubated for an additional 16 hours. Firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). The maximum induction fold was approximately [26.7]. Data are presented based on drug mass concentration.

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at 176 x g for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

FCGR2B(CD32B) [P31994-1](#)

MGILSFLPVLATESDWADCKSPQPWGHMLLWTA VLFLAPVAGTPAAPPKAVLKLEPQWINVLQEDSVTLTC
RGTHSPESDSIQWFHNGNLIPTHTQPSYRFKANNNDSGEYTCQTGQTSLSDPVHLTVLSEWLVLQTPHLEFQE
GETIVLRCHSWKDKPLVKVTFQNGKSKKFSRSDPNFNSIPQANHSHSGDYHCTGNIGYTLYSSKPVTITVQAPS
SSPMGIIVAVVTGIAVA AIVAAVVALIYCRKKRISALPGYPECREMGETLPEK PANPTNPDEADKVG AENTITY
SLLMHPDALEEPDDQNRI

Related Products

FcγR	
Cynomolgus_FcRn MDCK Cell Line	H_FCGR1A(CD64) CHO-K1 Cell Line

H_FCGR1A(CD64) HEK-293 Cell Line	H_FCGR2A(CD32A) CHO-K1 Cell Line
H_FCGR3A(CD16a) 158F CHO-K1 Cell Line	H_FCGR3A(CD16a) 158V CHO-K1 Cell Line
H_FCGR3B(CD16b) CHO-K1 Cell Line	H_FcRn CHO-K1 Cell Line
H_FcRn MDCK Cell Line	H_FcRn MDCK Cell Line
Mouse_FcRn MDCK Cell Line	
Anti-FcRn hIgG4 Reference Antibody(Rozabio)	Anti-H_FcRn IgG4 Antibody(Rozanolixizumab)
Anti-Mouse CD1632 mIgG2b Antibody(2.4G2)	
ADCCP	
ADCC FcγRIIIa(158F) Jurkat Effector Cell Line	ADCC FcγRIIIa(158V) DDX35TM Jurkat Effector Cell Line
ADCC FcγRIIIa(158V) Jurkat Effector Cell Line	ADCC M_FcγRIV Jurkat Effector Cell Line
ADCP FcγRIIa DDX35TM Jurkat Effector Cell Line	ADCP FcγRIIa Jurkat Effector Cell Line
ADCP FcγRIIa R131 Jurkat Effector Cell Line	ADCP FcγRIIb Jurkat Effector Cell Line
NK92	
ADCC FcγRIIIa(158V) NK-92 Effector Cell Line	

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.