

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

# **Product Sheet**

## Cynomolgus\_SLC34A2 HEK-293 Cell Line

Catalog number: GM-C21535

Version 3.3.1.241203

Cynomolgus\_SLC34A2 HEK-293 Cell Line is a clonal stable HEK-293 cell line that

**Description** constitutively expresses the cynomolgus SLC34A2 gene, constructed using lentiviral

technology.

**Quantity** 5E6 Cells per vial,1 mL

**Product Format** 1 vial of frozen cells

**Shipping** Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Cynomolgus\_SLC34A2

**Gene ID/Uniprot ID** XP\_005554673.1 (AA Met 1 - Leu 573)

Host Cell HEK-293

**Recovery Medium** DMEM+10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+0.75 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

**Growth Conditions** 37°C, 5% CO<sub>2</sub>

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

**Safety considerations** Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Materials**

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-SLC34A2 hIgG1 Antibody(XMT-1536)	Genomeditech/GM-26828AB

#### **Figures**

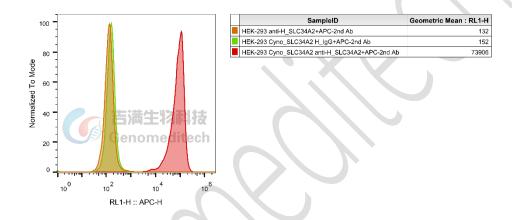


Figure 1 | Cynomolgus\_SLC34A2 HEK-293 Cell Line (Cat. GM-C21535) was determined by flow cytometry using Anti-H\_SLC34A2 hIgG1 Antibody(Cat. GM-26828AB).

### **Cell Recovery**

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

吉满生物科技 Genomeditech Genomeditech (Shanghai) Co.,Ltd.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

## **Cell Freezing**

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

## Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75  $\mu g/mL$  Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

#### **Notes**

a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.

b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

#### **Sequence**

SLC34A2 XP 005554673.1(ΔICD)

MAPWPELGDAQPNPSKYLEGATSQQPITPDKSKETNKNNTEVPVTKFELLPTYSTATLIEEPTEVDDPWNLPT LQDSGIKWSERDAKGKILCVFQGIGRLILLLGFLYFFVCSLDVLSSAFQLVGGKMAGQFFSNSSIMSNPLLGLV IGVLVTVFVQSSSTSTSIVVSMVASSLLTVRAAIPIIMGANIGTSITNTIVALMQVGDRSEFRRAFAGATVHDFF NWLSVLVLLPVEVATHYLEVVTQLIVESFHFKNGEDAPDLLKVITKPFTKLIVQLDKKVISQIAMNDETAKNK SLVKIWCKTFTNMTQMNVTVPSMANCTSPSLCWTDGIQIWTMKNVTYKENIAKCQHIFVNFHLPDLAVGIIL LIISLLVLCGCLIMIVKILGSVLKGQVATVIKKTINTDFPFPFAWLTGYLAILVGAGMTFIVQSSSVFTSALTPLI GIGVITIERAYPLTLGSNIGTTTTAILAALASPGNTLRSSLQIALCHFFFNISGILLWYPIPFTRLPIRMAKGLGNIS AKYRWFAVFYLIIFFFLIPLTVFGLSLAGWPVLVGVGVPVVFIIILVLCLRL

#### **Related Products**

SLC34A2(NaPi2b)		
Cynomolgus_SLC34A2 CHO-K1 Cell Line	H_SLC34A2 CHO-K1 Cell Line	
H_SLC34A2 HEK-293 Cell Line	H_SLC34A2 LLC1 Cell Line	
Mouse_SLC34A2 HEK-293 Cell Line	Rhesus_SLC34A2 CHO-K1 Cell Line	
Rhesus_SLC34A2 HEK-293 Cell Line		
Anti-SLC34A2 hIgG1 Antibody(XMT-1536)	Anti-SLC34A2 hIgG1 Reference Antibody (Upibio)	
$FOLR1(FR\alpha)$		
Cynomolgus_FOLR1 CHO-K1 Cell Line	H_FOLR1 CHO-K1 Cell Line	
H_FOLR1 CT26 Cell Line	H_FOLR1 HEK-293 Cell Line	
H_FOLR1 MC38 Cell Line		
Anti-FOLR1 hIgG1 Antibody(Mirvetuximab)	Anti-FOLR1 hIgG1 Reference Antibody(Mirvebio)	
Anti-FOLR1-DM4(Dar3.4)[Mirvetuximab soravtansine]		
Biotinylated Cynomolgus FOLR1 Protein; His-Avi Tag	Biotinylated Human FOLR1 Protein; His-Avi Tag	
Cynomolgus FOLR1 Protein; His Tag	Human FOLR1 Protein; His Tag	
CDH6		
Cynomolgus_CDH6 CHO-K1 Cell Line	H_CDH6 CHO-K1 Cell Line	
H_CDH6 HEK-293 Cell Line		
Anti-CDH6 hIgG1 Reference Antibody (Ralubio)	Anti-H_CDH6 hIgG1 Antibody(H01L02)	
ADC Related Product		
Anti-DXD Mouse IgG1 Antibody (23E21C5)	Anti-DXD Mouse IgG1 Antibody (4A5A12)	
Anti-Dxd Mouse IgG2a Antibody (17D6A4)	Anti-Eribulin Mouse IgG2a Antibody (10F8G4)	
Anti-MMAE Mouse IgG1 Antibody (11C10E3)	Anti-MMAE Mouse IgG2a Antibody (17A1K11)	
Anti-MMAE Mouse IgG2a Antibody (8F6A3)	Mouse anti Human IgG-MMAE(Dar4)	
Human IgG1 Isotype-DXD (Dar8)	Human IgG1 Isotype-Eribulin (Dar4)	
Human IgG1 Isotype-MMAE (Dar4)		
Recombinant DT3C Protein		



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

#### **Limited Use License Agreement**

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.