

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

Cynomolgus_BCMA HEK-293 Cell Line

Catalog number: GM-C24898

Version 3.3.1.251024

Cynomolgus_BCMA HEK-293 Cell Line is a clonal stable HEK-293 cell line that

Description constitutively expresses the cynomolgus BCMA gene, constructed using lentiviral

technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Cynomolgus_BCMA

Gene ID/Uniprot ID A0A2K5UD97

Host Cell HEK-293

Recovery Medium DMEM+10% FBS+1% P.S

Growth medium DMEM+10% FBS+1% P.S+0.75 μ g/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	ExCell/FSP500
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-BCMA hIgG1 Antibody(Belantamab)	Genomeditech/GM-52396AB

Figures

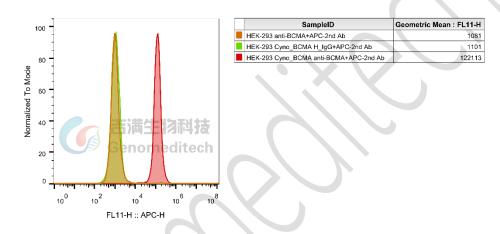


Figure 1 | Cynomolgus_BCMA HEK-293 Cell Line (Cat. GM-C24898) was determined by flow cytometry using Anti-BCMA hIgG1 Antibody(Belantamab) (Cat. GM-52396AB).

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- Aliquot 1 mL into each vial. c)
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor. c)
- Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell d) layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting. f)
- After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels. g)
- Incubate cultures at 37°C. h)

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Sequence

BCMA(TNFRSF17) A0A2K5UD97

 $\label{thm:construction} MLQMARQCSQNEYFDSLLHDCKPCQLRCSSTPPLTCQRYCNASMTNSVKGMNAILWTCLGLSLIISLAVFVL\\ TFLLRKMSSEPLKDEFKNTGSGLLGMANIDLEKGRTGDEIVLPRGLEYTVEECTCEDCIKNKPKVDSDHCFPL\\ PAMEEGATILVTTKTNDYCNSLSAALSVTEIEKSISAR$

Related Products

C	D40: CD40L	
H_CD40(TNFRSF5) Reporter 293 Cell Line	H_CD40(TNFRSF5) Reporter Jurkat Cell Line	
Cynomolgus_CD40 CHO-K1 Cell Line	Cynomolgus_CD40L CHO-K1 Cell Line	
H_CD40(TNFRSF5) CHO-K1 Cell Line	H_CD40(TNFRSF5) HEK-293 Cell Line	
H_CD40L CHO-K1 Cell Line	H_CD40L HEK-293 Cell Line	
Rabbit_CD40L NIH-3T3 Cell Line		
Anti-CD40L hIgG1 Reference Antibody (Frebio)	Anti-H_CD40 hIgG1 Antibody(APX005M)	
Anti-H_CD40 hIgG1 Antibody(ravagalimab)	Anti-H_CD40L hIgG1 Antibody(dapirolizumab)	
Anti-H_CD40L hIgG1 Antibody(frexalimab)		
Biotinylated Human CD40 Protein; His-Avi Tag	Cynomolgus CD40 Protein; His Tag	
Human CD40 Protein; His Tag	Human CD40L Protein; His Tag	
IFN-α		
IFNα Reporter HEK-293 Cell Line	IFNα Reporter MDCK Cell Line	
IFNα Reporter THP1 Cell Line		
BCM	IA:BAFFR:TACI	
H_BAFFR Jurkat Blockade Reporter Cell Line	H_BAFFR Reporter Cell Line	
H_BCMA Reporter Cell Line	H_TACI Reporter Cell Line	
Cynomolgus_BCMA CHO-K1 Cell Line	H_BCMA CHO-K1 Cell Line	
H_BCMA HEK-293 Cell Line		
Anti-BAFF hIgG1 Antibody(belimumab)	Anti-BAFFR hIgG1 Antibody(ianalumab)	
Anti-BCMA hIgG1 Antibody(Belantamab)	Anti-BCMA hIgG1 Antibody(SEA-BCMA)	
Anti-BCMA hIgG4 Antibody(BCMB69)		
Biotinylated Human BAFF Protein; His-Avi Tag	Cynomolgus BAFF Protein; His Tag	
Human BAFF Protein; His Tag	Mouse BAFF Protein; His Tag	
BD	CA2(CLEC4C)	
H_BDCA2 Reporter Jurkat Cell Line	Cynomolgus_BDCA2 CHO-K1 Cell Line	
Cynomolgus_BDCA2 Jurkat Cell Line	H_BDCA2 CHO-K1 Cell Line	
H_BDCA2 HEK-293 Cell Line	H_BDCA2 Jurkat Cell Line	
Anti-H_BDCA2 hIgG1 Antibody(Litifilimab)		
Cynomolgus BDCA2 Protein; His Tag	Human BDCA2 Protein; His Tag	
	CD3	
H_CD3D CD3E KO Jurkat Cell Line	ADCC FcyRIIIa(158V) Reporter Jurkat(CD3 KO) Cell Line	



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

CD3-CD2-tsAb Reporter Jurkat(CD58 KO) Cell Line	Jurkat CD3-BsAb Reporter Cell Line
Cynomolgus_CD3 HEK-293 Cell Line	Cynomolgus_CD3E(Membrane Bound ECD) CHO-K1 Cell Line
H_CD3 CHO-K1 Cell Line	H_CD3 HEK-293 Cell Line
H_CD3E(Membrane Bound ECD) CHO-K1 Cell Line	Mouse_CD3 HEK-293 Cell Line
Anti-CD3 epsilon hIgG1 Antibody [OKT-3 (muromonab)]	Anti-CD3 hIgG1 Antibody(CH2527)
Anti-mouse CD3ε mIgG2a Antibody(145-2C11)	

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly
 used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope;
 modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.