

Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

H_CD122 CD132 CHO-K1 Cell Line

Catalog number: GM-C26420

Version 3.3.1.250912

H_CD122 CD132 CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively **Description**

expresses the Human CD122 and CD132 genes, constructed using lentiviral technology.

Quantity 5E6 Cells per vial,1 mL

Product Format 1 vial of frozen cells

Shipping Shipped on dry ice

Storage Conditions Liquid nitrogen immediately upon receipt

Target Human_CD122 & Human_CD132

Gene ID/Uniprot ID P14784 & P31785-1

Host Cell CHO-K1

Recovery Medium F12K+10% FBS+1% P.S

Growth medium F12K+10% FBS+1% P.S+200 μg/mL G418+4 μg/mL Puromycin

Note None

Freezing Medium 90% FBS+10% DMSO

Growth properties Adherent

Growth Conditions 37°C, 5% CO₂

Mycoplasma Testing The cell line has been screened to confirm the absence of Mycoplasma species.

Safety considerations Biosafety Level 2

Note It is recommended to expand the cell culture and store a minimum of 10 vials at an early

passage for potential future use.




Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials


Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
G418	Genomeditech/GM-040402
Puromycin	Genomeditech/GM-040401
Anti-CD132(IL2RG) hIgG4 Antibody(REGN7257)	Genomeditech/GM-52334AB
Anti-CD122 hIgG1 Antibody(HuABC-2)	Genomeditech/GM-52319AB

Figures

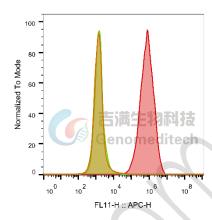

	SampleiD	Geometric Mean : FL11-H
	CHO-K1 anti-CD122+APC-2nd Ab	1407
	CHO-K1 H_CD122 CD132 H_lgG+APC-2nd Ab	1511
	CHO-K1 H_CD122 CD132 anti-CD122+APC-2nd Ab	7.19E5

Figure 1 | H_CD122 CD132 CHO-K1 Cell Line (Cat. GM-C26420) was determined by flow cytometry using Anti-CD122 hIgG1 Antibody(HuABC-2)(Cat. GM-52319AB)

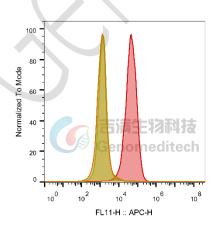

SampleID	Geometric Mean : FL11-H
CHO-K1 anti-CD132+APC-2nd Ab	1263
CHO-K1 H_CD122 CD132 H_lgG+APC-2nd Ab	1381
CHO-K1 H_CD122 CD132 anti-CD132+APC-2nd Ab	43232

Figure 2 | H_CD122 CD132 CHO-K1 Cell Line (Cat. GM-C26420) was determined by flow cytometry using Anti-CD132(IL2RG) hIgG4 Antibody(REGN7257)(Cat. GM-52334AB)



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+200 µg/mL G418+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- d) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.

 Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

IL2RB(CD122) P14784

MAAPALSWRLPLLILLLPLATSWASAAVNGTSQFTCFYNSRANISCVWSQDGALQDTSCQVHAWPDRRRWN QTCELLPVSQASWACNLILGAPDSQKLTTVDIVTLRVLCREGVRWRVMAIQDFKPFENLRLMAPISLQVVHV ETHRCNISWEISQASHYFERHLEFEARTLSPGHTWEEAPLLTLKQKQEWICLETLTPDTQYEFQVRVKPLQGE FTTWSPWSQPLAFRTKPAALGKDTIPWLGHLLVGLSGAFGFIILVYLLI*

IL2RG(CD132) P31785-1

MLKPSLPFTSLLFLQLPLLGVGLNTTILTPNGNEDTTADFFLTTMPTDSLSVSTLPLPEVQCFVFNVEYMNCTW
NSSSEPQPTNLTLHYWYKNSDNDKVQKCSHYLFSEEITSGCQLQKKEIHLYQTFVVQLQDPREPRRQATQML
KLQNLVIPWAPENLTLHKLSESQLELNWNNRFLNHCLEHLVQYRTDWDHSWTEQSVDYRHKFSLPSVDGQK
RYTFRVRSRFNPLCGSAQHWSEWSHPIHWGSNTSKENPFLFALEAVVISVGSMGLIISLLCVYFWLERTMPRIP
TLKNLEDLVTEYHGNFSAWSGVSKGLAESLQPDYSERLCLVSEIPPKGGALGEGPGASPCNQHSPYWAPPCY
TLKPET*

Related Products

Related Frouders				
IL-15				
H_IL15 Reporter Cell Line	Cynomolgus_CD122 HEK-293 Cell Line			
H_CD122 CHO-K1 Cell Line	H_CD122 HEK-293 Cell Line			
H_CD215(IL15RA) HEK-293 Cell Line				
IL-2				
H_CD122 CD132 Reporter Cell Line	H_CD25 CD122 CD132 Reporter Cell Line			
H_IL2 Reporter Cell Line	H_IL2 Reporter DDX35TM Cell Line			
Cynomolgus_CD25 HEK-293 Cell Line	H_CD25 CHO-K1 Cell Line			
H_CD25 HEK-293 Cell Line				
Anti-CD122 hIgG1 Antibody(HuABC-2)	Anti-CD132(IL2RG) hIgG4 Antibody(REGN7257)			
Anti-CD25 hIgG1 Antibody(Basiliximab)	Anti-mouse CD25 mIgG2a Antibody(PC-61.5.3)			
Anti-mouse CD25 RIgG1 Antibody(PC-61.5.3)				



Order: +86 021-68455258/50432826/50432825

Toll-free: +86 400 627 9288

Email: service@genomeditech.com

License Agreement:

By purchasing and using this cell line product, the user voluntarily agrees to accept and abide by the following policies:

- This cell line product is restricted to research use only and shall not be used for any commercial purposes.
- This product is strictly prohibited from being used in the diagnosis or treatment of human or animal diseases, and shall not be directly
 used in experiments involving humans.
- Users and their contractors engaged for their benefit may use this material and its derivatives only within the agreed research scope;
 modification of the material is not permitted, nor may it be distributed, sold, transferred, or otherwise provided to any other entity (including affiliates).
- If use beyond the above scope is required, prior written permission from Genomeditech (Shanghai) Co.,Ltd. must be obtained. For details, please contact Genomeditech (Shanghai) Co.,Ltd.