

# **Product Sheet**

## H\_CLDN18.2 MKN45 Cell Line(High Expression)

Catalog number: GM-C29870

Version 3.3.1.241120

Description	H_CLDN18.2 MKN45 Cell Line(High Expression) is a clonal stable MKN45 cell line that constitutively expresses the human CLDN18.2 gene, constructed using lentiviral technology.	
Quantity	5E6 Cells per vial,1 mL	
Product Format	1 vial of frozen cells	
Shipping	Shipped on dry ice	
Storage Conditions	Liquid nitrogen immediately upon receipt	
Target	Human_CLDN18.2 & C-3×Flag	
Gene ID/Uniprot ID	NP_001002026.1	
Host Cell	MKN45	
<b>Recovery Medium</b>	RPMI 1640+20% FBS+1% P.S	
Growth medium	RPMI 1640+20% FBS+1% P.S+1 µg/mL Puromycin	
Note	None	
Freezing Medium	90% FBS+10% DMSO	
Growth properties	Mixed: adherent and suspension	
Growth Conditions	37°C, 5% CO <sub>2</sub>	
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.	
Safety considerations	Biosafety Level 2	
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.	



Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Geometric Mean : FL11-H

3508


879

66943

#### **Materials**

Reagent	Manufacturer/Catalogue No.
RPMI 1640	VivaCell/C3010-0500
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-CLDN18.2 hIgG1 Antibody	In house/

#### Figures

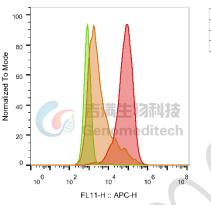



SampleID

MKN45 H\_CLDN18.2 anti-CLDN18.2+APC-2nd At

MKN45 anti-CLDN18.2+APC-2nd Ab

MKN45 H\_CLDN18.2 H\_IgG+APC-2nd Ab

Figure 1 | H\_CLDN18.2 MKN45 Cell Line(High Expression) (Cat. GM-C29870) was determined by flow cytometry using Anti-CLDN18.2 hIgG1 Antibody (In house).

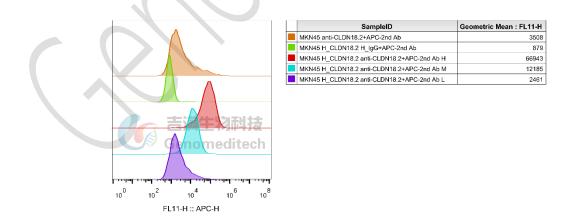


Figure 2 | H\_CLDN18.2 MKN45 Cell Line(High Expression) (Cat. GM-C29870) , H\_CLDN18.2 MKN45 Cell Line(Medium Expression) (Cat. GM-C29754) and H\_CLDN18.2 MKN45 Cell Line(Low Expression) (Cat. GM-C29874) were determined by flow cytometry using Anti-CLDN18.2 hIgG1 Antibody (In house).

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



#### **Cell Recovery**

Recovery Medium: RPMI 1640+20% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at  $-70^{\circ}$ C. Storage at  $-70^{\circ}$ C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO<sub>2</sub> in air atmosphere is recommended if using the medium described on this product sheet.

### **Cell Freezing**

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

#### Cell passage

Growth medium: RPMI 1640+20% FBS+1% P.S+1 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Under normal conditions, these cells exist as both adherent and round suspension cells.
- b) When changing the medium, take care to retain the suspension cells. During passaging, collect both the adherent and suspension cells together before subculturing.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 1 to 2 minutes at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.
  Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.



#### Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

#### Notes

- a) Under normal conditions, these cells exist as both adherent and round suspension cells.
- b) When changing the medium, take care to retain the suspension cells. During passaging, collect both the adherent and suspension cells together before subculturing.
- c) Once cell status stabilizes, the number of dead cells will decrease after passaging, the growth rate will become stable, cell morphology will be uniform, and the cells will appear robust.

#### Sequence

CLDN18.2(isoform2) NP\_001002026.1

MAVTACQGLGFVVSLIGIAGIIAATCMDQWSTQDLYNNPVTAVFNYQGLWRSCVRESSGFTECRGYFTLLGL PAMLQAVRALMIVGIVLGAIGLLVSIFALKCIRIGSMEDSAKANMTLTSGIMFIVSGLCAIAGVSVFANMLVTN FWMSTANMYTGMGGMVQTVQTRYTFGAALFVGWVAGGLTLIGGVMMCIACRGLAPEETNYKAVSYHAS GHSVAYKPGGFKASTGFGSNTKNKKIYDGGARTEDEVQSYPSKHDYVGSDYKDHDGDYKDHDIDYKDDDD K\*

#### **Related Products**

CLDN18				
Cynomolgus_CLDN18.2-eGFP CHO-K1 Cell Line	H_CLDN18(isoform2)-eGFP 293 Cell Line			
H_CLDN18.1-eGFP HEK-293 Cell Line	H_CLDN18.2 MC38 Cell Line			
H_CLDN18.2 MKN45 Cell Line	H_CLDN18.2 MKN45 Cell Line(Low Expression)			
H_CLDN18.2 MKN45 Cell Line(Medium Expression)	H_CLDN18.2(isoform2) CHO-K1 Cell Line			
H_CLDN18.2-eGFP CT-26 Cell Line	Mouse_CLDN18.2-eGFP CHO-K1 Cell Line			
Rat_CLDN18.2-eGFP CHO-K1 Cell Line	Rhesus_CLDN18.2-eGFP CHO-K1 Cell Line			
Anti-CLDN18.2 hIgG1 Antibody(LM-102)	Anti-CLDN18.2 hIgG1 Antibody(Zolbetuximab)			
HER3(ERBB3)				
Cynomolgus_ERBB3(HER3) CHO-K1 Cell Line	Cynomolgus_ERBB3(HER3) HEK-293 Cell Line			
H_ERBB3(HER3) CHO-K1 Cell Line	H_ERBB3(HER3) HEK-293 Cell Line			
H_ERBB3(HER3) MC38 Cell Line	Mouse_HER3(ERBB3) CHO-K1 Cell Line			
Anti-ERBB3(HER3) hIgG1 Reference Antibody(Patribio)	Anti-H_ERBB3(HER3) hIgG1 Antibody(Barecetamab)			
Human HER3 Protein; His Tag				
TROP2(TACSTD2)				
Cynomolgus_Trop2 CHO-K1 Cell Line	Cynomolgus_TROP2 HEK-293 Cell Line			
H_TROP2 CHO-K1 Cell Line	H_TROP2 CT26 Cell Line			
H_TROP2 HEK-293 Cell Line	H_TROP2 LLC1 Cell Line			
H_TROP2 MC38 Cell Line				

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!

Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

_	Email: service@genomeditech.com
Anti-H_TROP2 hIgG1 Antibody(Datopotamab)	Anti-TROP2 hIgG1 Antibody(Hu2G10-5)
Anti-Trop2 hIgG1 Reference Antibody (Sacbio)	Anti-Trop2 hIgG1 Reference Antibody(Datbio)
Anti-Trop2-DXD ADC(Dar4)[Datopotamab deruxtecan,Dato- DXD]	Anti-Trop2-SN38 ADC(Dar8)[Sacituzumab govitecan]
Human TROP2 Protein; His Tag	
GUCY2	C(GC-C)
H_GUCY2C CHO-K1 Cell Line	H_GUCY2C HEK-293 Cell Line
Anti-H_GUCY2C hIgG1 Antibody(Indusatumab)	
С	D3
Jurkat CD3-BsAb Reporter Cell Line	Cynomolgus_CD3 HEK-293 Cell Line
Cynomolgus_CD3E(Membrane Bound ECD) CHO-K1 Cell Line	H_CD3 CHO-K1 Cell Line
H_CD3 HEK-293 Cell Line	H_CD3E(Membrane Bound ECD) CHO-K1 Cell Line
Mouse_CD3 HEK-293 Cell Line	
Anti-CD3 epsilon hIgG1 Antibody [OKT-3 (muromonab)]	Anti-CD3 hIgG1 Antibody(CH2527)
CL	DN3
H_CLDN3 HEK-293 Cell Line	
Anti-CLDN3 hIgG1 Antibody(H4G3)	
CL	DN4
H_CLDN4 HEK-293 Cell Line	
Anti-CLDN4 hIgG1 Antibody(4B8)	
CL	DN6
Cynomolgus_CLDN6 CHO-K1 Cell Line	H_CLDN6 CHO-K1 Cell Line
H_CLDN6 HEK-293 Cell Line	H_CLDN6 LLC1 Cell Line
Mouse_CLDN6 CHO-K1 Cell Line	Rat_CLDN6 CHO-K1 Cell Line
Rhesus_CLDN6 CHO-K1 Cell Line	
Anti-Claudin6 hIgG1 Reference Antibody	Anti-CLDN6/9 hIgG1 Antibody
CL	DN9
Cynomolgus_CLDN9 CHO-K1 Cell Line	H_CLDN9 CHO-K1 Cell Line
H_CLDN9-eGFP HEK-293 Cell Line	
ADC Rela	ted Product
Anti-DXD Mouse IgG1 Antibody (23E21C5)	Anti-DXD Mouse IgG1 Antibody (4A5A12)
Anti-Dxd Mouse IgG2a Antibody (17D6A4)	Anti-Eribulin Mouse IgG2a Antibody (10F8G4)
Anti-MMAE Mouse IgG1 Antibody (11C10E3)	Anti-MMAE Mouse IgG2a Antibody (17A1K11)
Anti-MMAE Mouse IgG2a Antibody (8F6A3)	Mouse anti Human IgG-MMAE(Dar4)
Human IgG1 Isotype-DXD (Dar8)	Human IgG1 Isotype-Eribulin (Dar4)
Human IgG1 Isotype-DXD (Dar8) Human IgG1 Isotype-MMAE (Dar4)	Human IgG1 Isotype-Eribulin (Dar4)



#### **Limited Use License Agreement**

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

吉満生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd 上海市浦东新区康威路 299 号1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!