

Product Sheet

H_CEACAM5(CEA) CT26 Cell Line

Catalog number: GM-C31428

Version 3.3.1.241112

Description	H_CEACAM5(CEA) CT26 Cell Line is a clonal stable CT26 cell line constitutively expressing human CEACAM5(CEA).	
Quantity	5E6 Cells per vial,1 mL	
Product Format	1 vial of frozen cells	
Shipping	Shipped on dry ice	
Storage Conditions	Liquid nitrogen immediately upon receipt	
Target	Human_CEACAM5(CEA)	
Gene ID/Uniprot ID	P06731-1	
Host Cell	CT26	
Recovery Medium	RPMI 1640+10% FBS+1% P.S	
Growth medium	RPMI 1640+10% FBS+1% P.S+4 µg/mL Blasticidin+4 µg/mL Puromycin	
Note	None	
Freezing Medium	90% FBS+10% DMSO	
Growth properties	Adherent	
Growth Conditions	37°C, 5% CO ₂	
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.	
Safety considerations	Biosafety Level 2	
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.	



Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
RPMI 1640	gibco/C11875500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/GM-040404
Puromycin	Genomeditech/GM-040401
Anti-CEA hIgG1 Antibody(Cibisatamab)	Genomeditech/GM-28531AB

Figures

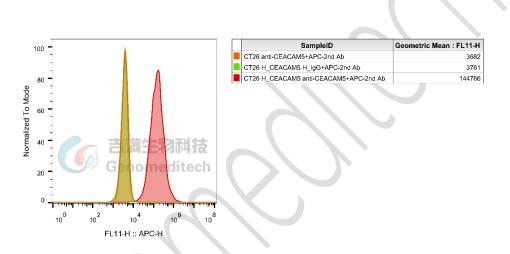


Figure 1 | H_CEACAM5(CEA) CT26 Cell Line (Cat. GM-C31428) was determined by flow cytometry using Anti-CEA hIgG1 Antibody(Cibisatamab) (Cat. GM-28531AB).

Cell Recovery

Recovery Medium: RPMI 1640+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70° C. Storage at -70° C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: RPMI 1640+10% FBS+1% P.S+4 µg/mL Blasticidin+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.
 Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Sequence

CEACAM5(CEA) P06731-1

MESPSAPPHRWCIPWQRLLLTASLLTFWNPPTTAKLTIESTPFNVAEGKEVLLLVHNLPQHLFGYSWYKGER VDGNRQIIGYVIGTQQATPGPAYSGREIIYPNASLLIQNIIQNDTGFYTLHVIKSDLVNEEATGQFRVYPELPKP SISSNNSKPVEDKDAVAFTCEPETQDATYLWWVNNQSLPVSPRLQLSNGNRTLTLFNVTRNDTASYKCETQN



PVSARRSDSVILNVLYGPDAPTISPLNTSYRSGENLNLSCHAASNPPAQYSWFVNGTFQQSTQELFIPNITVNN SGSYTCQAHNSDTGLNRTTVTTITVYAEPPKPFITSNNSNPVEDEDAVALTCEPEIQNTTYLWWVNNQSLPVS PRLQLSNDNRTLTLLSVTRNDVGPYECGIQNELSVDHSDPVILNVLYGPDDPTISPSYTYYRPGVNLSLSCHAA SNPPAQYSWLIDGNIQQHTQELFISNITEKNSGLYTCQANNSASGHSRTTVKTITVSAELPKPSISSNNSKPVED KDAVAFTCEPEAQNTTYLWWVNGQSLPVSPRLQLSNGNRTLTLFNVTRNDARAYVCGIQNSVSANRSDPVT LDVLYGPDTPIISPPDSSYLSGANLNLSCHSASNPSPQYSWRINGIPQQHTQVLFIAKITPNNNGTYACFVSNLA TGRNNSIVKSITVSASGTSPGLSAGATVGIMIGVLVGVALI

Related Products

HAVCR2(T	IM3): CEACAM1
Cynomolgus_HAVCR2(TIM3) CHO-K1 Cell Line	H_CEACAM1(CD66a) CHO-K1 Cell Line
H_HAVCR2(TIM3) CHO-K1 Cell Line	H_HAVCR2(TIM3) HEK-293 Cell Line
Anti-H_CEACAM1(CD66a) hIgG4 Antibody(MK-6018)	Anti-H_HAVCR2(TIM3) hIgG4 Antibody(Sabatolimab)
Anti-TIM3 hIgG2 Antibody(M6903)	
CH	EACAM5
Cynomolgus_CEACAM5(CEA) CHO-K1 Cell Line	Cynomolgus_CEACAM5(CEA) HEK-293 Cell Line
H_CEACAM5(CEA) CHO-K1 Cell Line	H_CEACAM5(CEA) HEK-293 Cell Line
H_CEACAM5(CEA) MC38 Cell Line	
Anti-CEA hIgG1 Antibody(Cibisatamab)	Anti-CEACAM5 hIgG1 Reference Antibody (Labbio)
Anti-CEACAM5 hIgG1 Reference Antibody (Tusbio)	
Human CEACAM5 Protein; His Tag	
CH	EACAM8
Anti-CEACAM8 hIgG1 Antibody(Besilesomab)	
CH	EACAM6
Cynomolgus_CEACAM6(CD66c) CHO-K1 Cell Line	H_CEACAM6(CD66c) CHO-K1 Cell Line
H_CEACAM6(CD66c) HEK-293 Cell Line	H_CEACAM6(CD66c) HEK-293 Cell Line (High Expression)
Anti-H_CEACAM6(CD66c) hIgG1 Antibody(Tinurilimab)	
CH	EACAM7
Cynomolgus_CEACAM7 CHO-K1 Cell Line	H_CEACAM7 CHO-K1 Cell Line
H_CEACAM7 HEK-293 Cell Line	
Anti-CEACAM7 hIgG4 Antibody(2869)	
ADC R	elated Product
Anti-DXD Mouse IgG1 Antibody (23E21C5)	Anti-DXD Mouse IgG1 Antibody (4A5A12)
Anti-Dxd Mouse IgG2a Antibody (17D6A4)	Anti-Eribulin Mouse IgG2a Antibody (10F8G4)
Anti-MMAE Mouse IgG1 Antibody (11C10E3)	Anti-MMAE Mouse IgG2a Antibody (17A1K11)
Anti-MMAE Mouse IgG2a Antibody (8F6A3)	Mouse anti Human IgG-MMAE(Dar4)
Human IgG1 Isotype-DXD (Dar8)	Human IgG1 Isotype-Eribulin (Dar4)
Human IgG1 Isotype-MMAE (Dar4)	
Recombinant DT3C Protein	

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd 上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!