

Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

Mouse_CD3 HEK-293 Cell Line

Catalog number: GM-C33386

Version 3.3.1.241017

The CD3 complex, also known as the T3 complex, is a multimeric protein complex composed of four different polypeptide chains: epsilon (ϵ), gamma (γ), delta (δ), and zeta (ζ). These chains assemble and function as three pairs of dimers: $\epsilon \gamma$, $\epsilon \delta$, and $\zeta \zeta$. The CD3 protein complex is a defining feature of the T-cell lineage, thus anti-CD3 antibodies can be effectively used as markers for T cells.

Mouse_CD3 HEK-293 Cell Line is a clonal stable HEK-293 cell line constitutively expressing mouse CD3 complex.





Specifications

Quantity	5E6 Cells per vial,1 mL	
Product Format	1 vial of frozen cells	
Shipping	Shipped on dry ice	
Storage Conditions	Liquid nitrogen immediately upon receipt	
Recovery Medium	DMEM+10% FBS+1% P.S	
Growth medium	DMEM+10% FBS+1% P.S+400 µg/mL G418+125 µg/mL Hygromycin	
Note	None	
Freezing Medium	90% FBS+10% DMSO	
Growth properties	Adherent	
Growth Conditions	37°C, 5% CO ₂	
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.	
Safety considerations	Biosafety Level 2	
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.	

Materials

Reagent	Manufacturer/Catalogue No.
Hygromycin	Genomeditech/GM-040403
G418	Genomeditech/GM-040402
Pen/Strep	Thermo/15140-122
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
DMEM	Gibco/C11995500BT
Hieff® qPCR SYBR Green Master Mix (Low Rox Plus)	Yeasen/11202ES08
Hifair® II 1st Strand cDNA Synthesis SuperMix	Yeasen/11120ES60



Figures

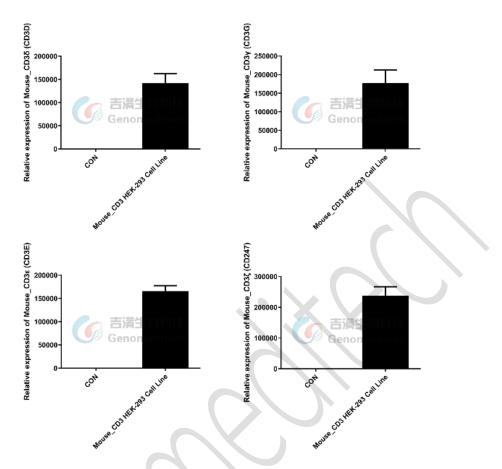


Figure 1 | The mRNA expression levels of Mouse_CD3D, CD3G, CD3E and CD247 in the Mouse_CD3 HEK-293 Cell Line (Cat. GM-C33386) were determined by RT-Qpcr.

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70° C. Storage at -70° C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+400 µg/mL G418+125 µg/mL Hygromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.
 Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.

Related Products

CD28			
H_CD28 Reporter Jurkat Cell Line	Cynomolgus_CD28 CHO-K1 Cell Line		
H_CD28 CHO-K1 Cell Line	H_CD28 HEK-293 Cell Line		
Anti-CD28 hIgG4 Antibody(FR104)	Anti-H_CD28 hIgG4 Antibody(Theralizumab)		
CD3			
Jurkat CD3-BsAb Reporter Cell Line	Cynomolgus_CD3 HEK-293 Cell Line		

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Cynomolgus_CD3E(Membrane Bound ECD) CHO-K1 Cell Line	H_CD3 CHO-K1 Cell Line
H_CD3 HEK-293 Cell Line	H_CD3E(Membrane Bound ECD) CHO-K1 Cell Line
Anti-CD3 epsilon hIgG1 Antibody [OKT-3 (muromonab)]	Anti-CD3 hIgG1 Antibody(CH2527)

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.