

# Product Sheet

## H\_CD40L CHO-K1 Cell Line

Catalog number: GM-C35532

Version 3.3.1.241128

<b>Description</b>	H_CD40L CHO-K1 Cell Line is a clonal stable CHO-K1 cell line that constitutively expresses the human CD40 gene, constructed using lentiviral technology.
<b>Quantity</b>	5E6 Cells per vial, 1 mL
<b>Product Format</b>	1 vial of frozen cells
<b>Shipping</b>	Shipped on dry ice
<b>Storage Conditions</b>	Liquid nitrogen immediately upon receipt
<b>Target</b>	Human_CD40L
<b>Gene ID/Uniprot ID</b>	P29965
<b>Host Cell</b>	CHO-K1
<b>Recovery Medium</b>	F12K+10% FBS+1% P.S
<b>Growth medium</b>	F12K+10% FBS+1% P.S+4 µg/mL Puromycin
<b>Note</b>	None
<b>Freezing Medium</b>	90% FBS+10% DMSO
<b>Growth properties</b>	Adherent
<b>Growth Conditions</b>	37°C, 5% CO <sub>2</sub>
<b>Mycoplasma Testing</b>	The cell line has been screened to confirm the absence of Mycoplasma species.
<b>Safety considerations</b>	Biosafety Level 2
<b>Note</b>	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

## Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-H_CD40L hIgG1 Antibody(dapirolizumab)	Genomeditech/GM-83977AB

## Figures

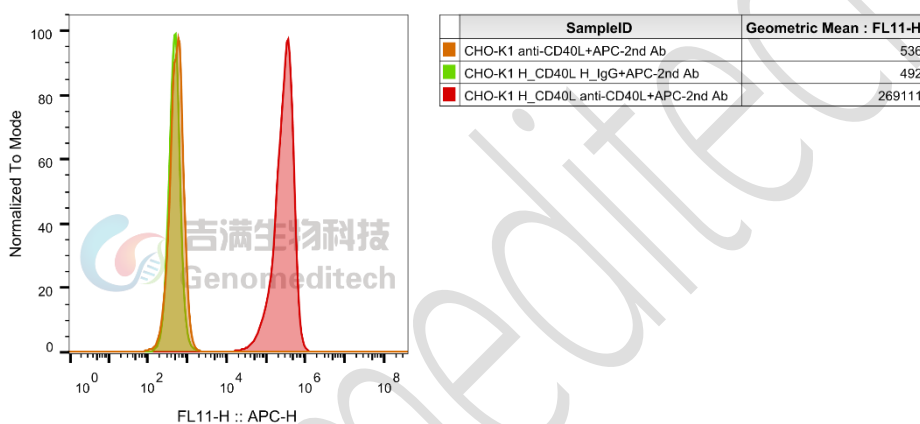


Figure 1 | H\_CD40L CHO-K1 Cell Line (Cat. GM-C35532) was determined by flow cytometry using Anti-H\_CD40L hIgG1 Antibody(dapirolizum) (Cat. [GM-83977AB](#)).

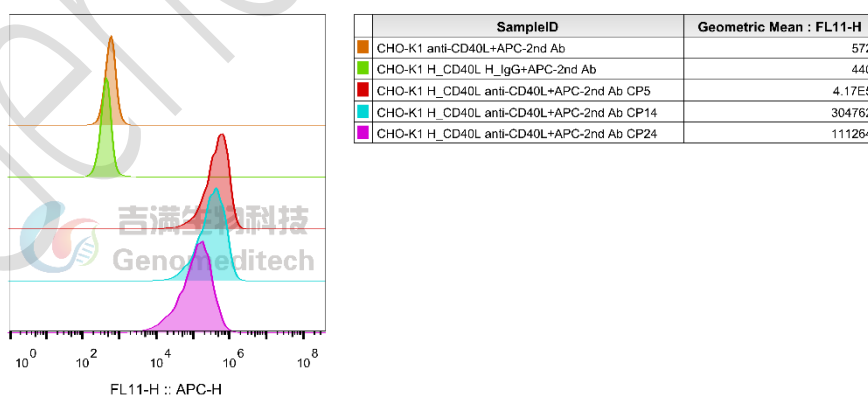


Figure 2 | The passage stability of the H\_CD40L CHO-K1 Cell Line (Cat. GM-C35532) was determined by flow cytometry using Anti-H\_CD40L hIgG1 Antibody(dapirolizum) (Cat. [GM-83977AB](#)).

## Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at  $-70^{\circ}\text{C}$ . Storage at  $-70^{\circ}\text{C}$  will result in loss of viability.

- Thaw the vial by gentle agitation in a  $37^{\circ}\text{C}$  water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately  $176 \times g$  for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- Incubate the culture at  $37^{\circ}\text{C}$  in a suitable incubator. A 5%  $\text{CO}_2$  in air atmosphere is recommended if using the medium described on this product sheet.

## Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at  $176 \times g$  for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to  $5 \times 10^6$  cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at  $-80^{\circ}\text{C}$  for at least 1 day, then transfer to liquid nitrogen as soon as possible.

## Cell passage

Growth medium: F12K+10% FBS+1% P.S+4  $\mu\text{g}/\text{mL}$  Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- Remove and discard culture medium.
- Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at  $37^{\circ}\text{C}$ ).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at  $37^{\circ}\text{C}$  to facilitate dispersal.
- Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- Incubate cultures at  $37^{\circ}\text{C}$ .

**Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended**

---

**Medium Renewal: Every 2 to 3 days**

## Notes

- a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

## Sequence

CD40L(CD40LG) [P29965](#)

MIETYNQTSRPSAATGLPISMKIFMYLLTVFLITQMIGSALFAVYLHRRLDKIEDERNLHEDFVFMKTIQRCNT  
 GERSLSLLNCEEIKSQFEGFVKDIMLNKEETKKENSFEMQKGDQNPQIAAHVISEASSKTTSVLQWAEKGYYT  
 MSNNLVTLENGKQLTVKRQGLYYIYAQVTFCSNREASSQAPFIASLCLKSPGRFERILLRAANTHSSAKPCGQ  
 QSIHLGGVFEHQPGASVFNVTDPSSQVSHGTGFTSFGLLKL

## Related Products

CD40: CD40L	
<a href="#">H_CD40(TNFRSF5) Reporter 293 Cell Line</a>	<a href="#">H_CD40(TNFRSF5) Reporter Jurkat Cell Line</a>
<a href="#">Cynomolgus_CD40 CHO-K1 Cell Line</a>	<a href="#">Cynomolgus_CD40L CHO-K1 Cell Line</a>
<a href="#">H_CD40(TNFRSF5) CHO-K1 Cell Line</a>	<a href="#">H_CD40(TNFRSF5) HEK-293 Cell Line</a>
<a href="#">H_CD40L HEK-293 Cell Line</a>	
<a href="#">Anti-H_CD40 hIgG1 Antibody(APX005M)</a>	<a href="#">Anti-H_CD40 hIgG1 Antibody(ravagalimab)</a>
<a href="#">Anti-H_CD40L hIgG1 Antibody(dapirolizumab)</a>	<a href="#">Anti-H_CD40L hIgG1 Antibody(frexalimab)</a>
<a href="#">Biotinylated Human CD40 Protein; His-Avi Tag</a>	<a href="#">Cynomolgus CD40 Protein; His Tag</a>
<a href="#">Human CD40 Protein; His Tag</a>	<a href="#">Human CD40L Protein; His Tag</a>
IFN- $\alpha$	
<a href="#">IFN<math>\alpha</math> Reporter HEK-293 Cell Line</a>	<a href="#">IFN<math>\alpha</math> Reporter MDCK Cell Line</a>
<a href="#">IFN<math>\alpha</math> Reporter THP1 Cell Line</a>	
BCMA:BAFFR:TACI	
<a href="#">H_BAFFR Reporter Cell Line</a>	<a href="#">H_BCMA Reporter Cell Line</a>
<a href="#">H_TACI Reporter Cell Line</a>	<a href="#">Cynomolgus_BCMA CHO-K1 Cell Line</a>
<a href="#">H_BCMA CHO-K1 Cell Line</a>	<a href="#">H_BCMA HEK-293 Cell Line</a>
<a href="#">Anti-BAFF hIgG1 Antibody(belumumab)</a>	<a href="#">Anti-BAFFR hIgG1 Antibody(ianalumab)</a>
<a href="#">Anti-BCMA hIgG1 Antibody(Belantamab)</a>	<a href="#">Anti-BCMA hIgG1 Antibody(SEA-BCMA)</a>
<a href="#">Anti-BCMA hIgG4 Antibody(BCMB69)</a>	
<a href="#">Biotinylated Human BAFF Protein; His-Avi Tag</a>	<a href="#">Cynomolgus BAFF Protein; His Tag</a>
<a href="#">Human BAFF Protein; His Tag</a>	<a href="#">Mouse BAFF Protein; His Tag</a>
BDCA2(CLEC4C)	
<a href="#">H_BDCA2 Reporter Jurkat Cell Line</a>	<a href="#">Cynomolgus_BDCA2 CHO-K1 Cell Line</a>
<a href="#">Cynomolgus_BDCA2 Jurkat Cell Line</a>	<a href="#">H_BDCA2 CHO-K1 Cell Line</a>
<a href="#">H_BDCA2 HEK-293 Cell Line</a>	<a href="#">H_BDCA2 Jurkat Cell Line</a>

Anti-H_BDCA2 hIgG1 Antibody(Litifilimab)	
Cynomolgus BDCA2 Protein; His Tag	Human BDCA2 Protein; His Tag

## Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

Genomeditech