

Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Product Sheet

H_MUSK HEK-293 Cell Line

Catalog number: GM-C37364

Version 3.3.1.241212

Description	H_MUSK HEK-293 Cell Line is a clonal stable HEK-293 cell line that constitutively expresses the human MUSK gene, constructed using lentiviral technology.	
Quantity	5E6 Cells per vail, 1 mL	
Product Format	1 vial of frozen cells	
Shipping	Shipped on dry ice	
Storage Conditions	Liquid nitrogen immediately upon receipt	
Target	Human MUSK	
Gene ID/Uniprot ID	O15146-1	
Host Cell	HEK-293	
Recovery Medium	DMEM+10% FBS+1% P.S	
Growth medium	DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin	
Note	None	
Freezing Medium	90% FBS+10% DMSO	
Growth properties	Adherent	
Growth Conditions	37°C, 5% CO ₂	
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.	
Safety considerations	Biosafety Level 2	
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.	



Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Puromycin	Genomeditech/GM-040401
Anti-MUSK hIgG1 Antibody(3B2g2m1)	Genomeditech/GM-87838AB

Figures

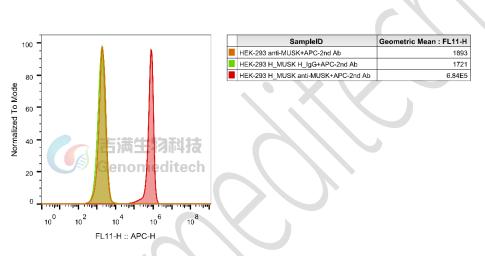


Figure 1 | H_MUSK HEK-293 Cell Line (Cat. GM-C37364) was determined by flow cytometry using Anti-MUSK hIgG1 Antibody(3B2g2m1) (Cat. GM-87838AB)

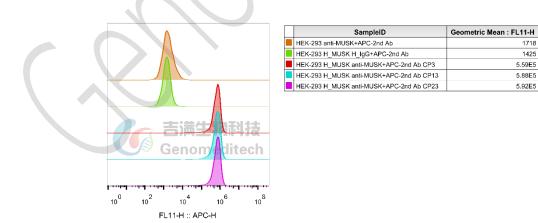


Figure 2 | The passage stability of the H_MUSK HEK-293 Cell Line (Cat. GM-C37364) was determined by flow cytometry using Anti-MUSK hIgG1 Antibody(3B2g2m1) (Cat. GM-87838AB).

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1, No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70° C. Storage at -70° C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.
- d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+0.75 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.
 Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途, 严禁用于人体治疗! For research use only!



- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

Sequence

MUSK 015146-1

MRELVNIPLVHILTLVAFSGTEKLPKAPVITTPLETVDALVEEVATFMCAVESYPQPEISWTRNKILIKLFDTRY SIRENGQLLTILSVEDSDDGIYCCTANNGVGGAVESCGALQVKMKPKITRPPINVKIIEGLKAVLPCTTMGNPK PSVSWIKGDSPLRENSRIAVLESGSLRIHNVQKEDAGQYRCVAKNSLGTAYSKVVKLEVEVFARILRAPESHN VTFGSFVTLHCTATGIPVPTITWIENGNAVSSGSIQESVKDRVIDSRLQLFITKPGLYTCIATNKHGEKFSTAKA AATISIAEWSKPQKDNKGYCAQYRGEVCNAVLAKDALVFLNTSYADPEEAQELLVHTAWNELKVVSPVCRP AAEALLCNHIFQECSPGVVPTPIPICREYCLAVKELFCAKEWLVMEEKTHRGLYRSEMHLLSVPECSKLPSMH WDPTACARLPHLDYNKENLKTFPPMTSSKPSVDIPNLPSSSSSSFSVSPTYSMTVIISIMSSFAIFVLLTITTLYCC RRRKQWKNKKRESAAVTLTTLPSELLLDRLHPNPMYQRMPLLLNPKLLSLEYPRNNIEYVRDIGEGAFGRVF QARAPGLLPYEPFTMVAVKMLKEEASADMQADFQREAALMAEFDNPNIVKLLGVCAVGKPMCLLFEYMAY GDLNEFLRSMSPHTVCSLSHSDLSMRAQVSSPGPPPLSCAEQLCIARQVAAGMAYLSERKFVHRDLATRNCL VGENMVVKIADFGLSRNIYSADYYKANENDAIPIRWMPPESIFYNRYTTESDVWAYGVVLWEIFSYGLQPYY GMAHEEVIYYVRDGNILSCPENCPVELYNLMRLCWSKLPADRPSFTSIHRILERMCERAEGTVSV

Related Products

MUSK		
Cynomolgus_MUSK HEK-293 Cell Line	H_MUSK CHO-K1 Cell Line	
Mouse_MUSK HEK-293 Cell Line	Rat_MUSK HEK-293 Cell Line	
Anti-MUSK hIgG1 Antibody(3B2g2m1)		
Anti-MUSK mIgG2a Antibody(3B2g2m1)		
Cynomolgus MUSK Protein; hFc Tag	Cynomolgus MUSK Protein; His Tag	
Human MUSK Protein; hFc Tag	Human MUSK Protein; His Tag	
Mouse MUSK Protein; hFc Tag	Mouse MUSK Protein; His Tag	

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd 上海市浦东新区康威路 299 号1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

吉満生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd 上海市浦东新区康威路 299 号1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!