

Product Sheet

H_CALCR Reporter CHO-K1 Cell Line


Catalog number: GM-C37605

Version 3.3.1.250124

CALCR(Calcitonin Receptor) is a G protein-coupled receptor encoded by the CALCR gene, mainly found in bones, kidneys, and the central nervous system. It regulates calcium metabolism and bone homeostasis by inhibiting osteoclast activity and promoting calcium excretion, playing a significant role in conditions like osteoporosis.

When calcitonin binds to CALCR, it activates different signaling pathways through G protein coupling mechanisms. Activates Gs protein, increasing adenylate cyclase activity and cAMP levels, which inhibits osteoclasts and reduces bone resorption. Stimulates phospholipase C to release intracellular calcium, enhancing physiological responses and regulating bone metabolism and inflammation. In summary, CALCR is vital for calcium homeostasis and may be a therapeutic target for related diseases.

H_CALCR Reporter CHO-K1 Cell Line is a clonal stable cell line constructed using lentiviral technology, constitutive expression of the CALCR gene, along with signal-dependent expression of a luciferase reporter gene. When calcitonin binds to CALCR, it activates downstream signaling pathways, leading to the expression of luciferase. The luciferase activity measurement indicates the activation level of the signaling pathway and can thus be used to evaluate the in vitro effects of drugs related to CALCR.

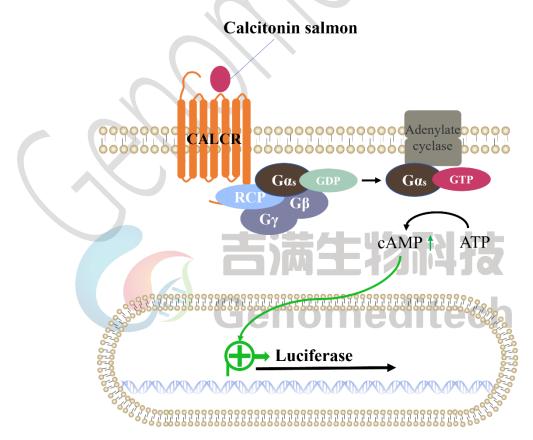

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途, 严禁用于人体治疗! For research use only!



Specifications

Quantity	5E6 Cells per vial,1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt
Recovery Medium	F12K+10% FBS+1% P.S
Growth medium	F12K+10% FBS+1% P.S+4 µg/mL Blasticidin+4 µg/mL Puromycin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂
Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
F12K	BOSTER/PYG0036
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/GM-040404
Puromycin	Genomeditech/GM-040401
Calcitonin salmon (Salmon calcitonin)	GlpBio/GC32851
GMOne-Step Luciferase Reporter Gene Assay Kit	Genomeditech/GM-040503



Figures

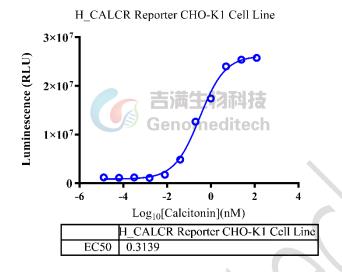


Figure 1 | Response to Calcitonin salmon. The H_CALCR Reporter CHO-K1 Cell Line (Cat. GM-C37605) at a concentration of 1E4 cells/well (96-well format) was stimulated with serial dilutions of Calcitonin salmon (GLPBIO/GC32851) in assay buffer (F12K + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). The maximum induction fold was approximately [22.4]. Data are shown by drug molar concentration.

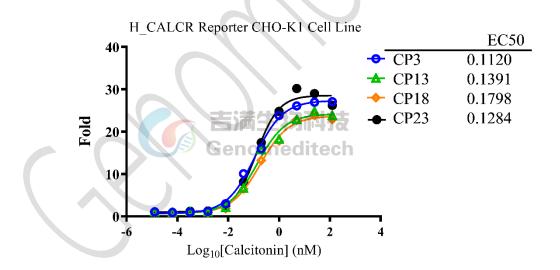


Figure 2 | The passage stability of response to Calcitonin salmon. The passage 3, 13, 18, and 23 of H_CALCR Reporter CHO-K1 Cell Line (Cat. GM-C37605) at a concentration of 1E4 cells/well (96-well format) were stimulated with serial dilutions of Calcitonin salmon (GLPBIO/GC32851) in assay buffer (F12K + 1% FBS + 1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. GM-040503). Data are shown by drug molar concentration.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!



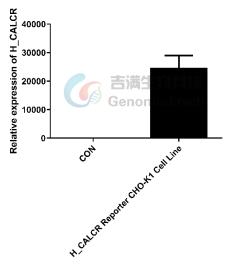


Figure 3 | The mRNA expression levels of H_CALCR in the H_CALCR Reporter CHO-K1 Cell Line (Cat. GM-C37605) were determined by RT-qPCR.

Cell Recovery

Recovery Medium: F12K+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70° C. Storage at -70° C will result in loss of viability.

- a) Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 3 minutes).
- b) Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- c) Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- d) Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- e) Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- a) Centrifuge at 176 x g for 3 minutes to collect cells.
- b) Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- c) Aliquot 1 mL into each vial.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only !



d) Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: F12K+10% FBS+1% P.S+4 µg/mL Blasticidin+4 µg/mL Puromycin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Remove and discard culture medium.
- b) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- c) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 2 to 3 minutes at 37°C).
- Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach.
 Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- e) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- f) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- g) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:4 - 1:5 is recommended

Medium Renewal: Every 2 to 3 days

Notes

a) After the stabilization of the cell condition, there will be fewer dead cells post-passage, the cell growth rate will tend to stabilize, cell morphology will become uniform, and the cells will appear robust.

Related Products

GCGR			
H_GCGR Reporter CHO-K1 Cell Line	H_GCGR Reporter HEK-293 Cell Line		
H_GCGR Reporter HEK-293 DDX35TM Cell Line	H_GCGR CHO-K1 Cell Line		
H_GCGR HEK-293 Cell Line	Mouse_GCGR HEK-293 Cell Line		
Anti-H_GCGR hIgG2 Antibody(volagidemab)			
GLP1R			
H_GLP1R Reporter CHO-K1 Cell Line	H_GLP1R Reporter HEK-293 Cell Line		
H_GLP1R Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GLP1R HEK-293 Cell Line		
H_GLP1R CHO-K1 Cell Line	H_GLP1R HEK-293 Cell Line		
Mouse_GLP1R HEK-293 Cell Line			
Anti-GLP1R hIgG1 Antibody(mAb-36986)	Anti-H_GLP1R hIgG1 Antibody(glutazumab)		
FGF21:FGFR			
H_FGF21 Reporter HEK-293 Cell Line			
CALCA(CGRP): CALCRL RAMP			

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1,No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only!

Genomeditech

Genomeditech (Shanghai) Co.,Ltd. Order: +86 021-68455258/50432826/50432825 Toll-free: +86 400 627 9288 Email: service@genomeditech.com

H_CALCRL RAMP1 Reporter HEK-293 Cell Line	H_CALCRL RAMP1 Reporter HEK-293 DDX35TM Cell Line			
Cynomolgus_CALCRL RAMP1 HEK-293 Cell Line	H_CALCRL RAMP1 CHO-K1 Cell Line			
H_CALCRL RAMP1 HEK-293 Cell Line				
Anti-CALCRL RAMP1 hIgG2 Antibody(Erenumab)				
GIP:GIPR				
H_GIPR Reporter CHO-K1 Cell Line	H_GIPR Reporter HEK-293 Cell Line			
H_GIPR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GIPR HEK-293 Cell Line			
H_GIPR CHO-K1 Cell Line	H_GIPR HEK-293 Cell Line			
Mouse_GIPR HEK-293 Cell Line				
Anti-H_GIPR hIgG1 Antibody(AMG-133)				
ACVR2A: A	ACTRIIB: Active A			
ACVR2A KO HEK-293 Cell Line	Activin A Reporter Cell Line			
H_ACVR2A Reporter Cell Line	H_ACVR2B Reporter Cell Line			
ACVR2B KO HEK-293 Cell Line	H_ACVR2A HEK-293(ACVR2B KO) Cell Line			
H_ACVR2B CHO-K1 Cell Line	H_ACVR2B HEK-293(ACVR2A KO) Cell Line			
Anti-ACVR2B hIgG1 Antibody(Bimagrumab)	Anti-ACVR2B hIgG1 Antibody(Fab-17G05)			
Anti-ACVR2B mIgG2a Antibody(Bimagrumab)	Anti-H_ACVR2B hIgG1 Reference Antibody(Bimbio)			
Biotinylated Human ACVR2A Protein; His-Avi Tag	Biotinylated Human ACVR2B Protein; His-Avi Tag			
Biotinylated Mouse ACVR2A Protein; His-Avi Tag	Biotinylated Mouse ACVR2B Protein; His-Avi Tag			
Human Activin A Protein; His Tag	Human Activin B Protein; His Tag			
Human ACVR2A Protein; hFc Tag	Human ACVR2A Protein; His Tag			
Human ACVR2B Protein; hFc Tag	Human ACVR2B Protein; His Tag			
Human latent GDF-8 Protein; His Tag	Mouse ACVR2B Protein; His Tag			
AMY: CALCR RAMP				
H_CALCR RAMP3(AMY3) Reporter CHO-K1 Cell Line				

H_CALCR RAMP3(AMY3) Reporter CHO-K1 Cell Line

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

吉满生物科技(上海)有限公司 Genomeditech (Shanghai) Co., Ltd

上海市浦东新区康威路 299 号 1 幢东区 505-507 邮编 201315 505-507,5th Floor, East District, Building 1, No.299 Kangwei Road, Pudong New Area, Shanghai 本公司产品仅供科研用途,严禁用于人体治疗! For research use only !