

Product Sheet


BRE Reporter 293 Cell Line

Catalog number: GM-C37665

Version 3.3.1.250227

Bone Morphogenetic Proteins (BMPs) are important members of the Transforming Growth Factor- β (TGF- β) superfamily, initially discovered for their role in bone and cartilage formation. The BMP signaling pathway plays a critical role in various biological processes, including embryonic development, cell differentiation, proliferation, apoptosis, as well as tissue regeneration and repair.

BRE Reporter 293 Cell Line is a clonal stable cell line with signal-dependent expression of a luciferase reporter gene stably transfected with the BMP Response Element (BRE) using lentiviral technology, and it endogenously expresses Type I Reporter and Type II Reporter gene. When ligands such as BMP-2 bind to their receptors, they activate downstream signaling pathways, which subsequently drive the expression of luciferase. By measuring the fluorescence signal, the expression of luciferase can be determined, indirectly evaluating the effect of the ligand on the signaling pathway. The luciferase readout represents the activation level of the signaling pathway. Therefore, this cell line is widely used for the screening of BMP ligands, receptor agonists, or inhibitors, as well as for studying the regulatory mechanisms of the signaling pathway.

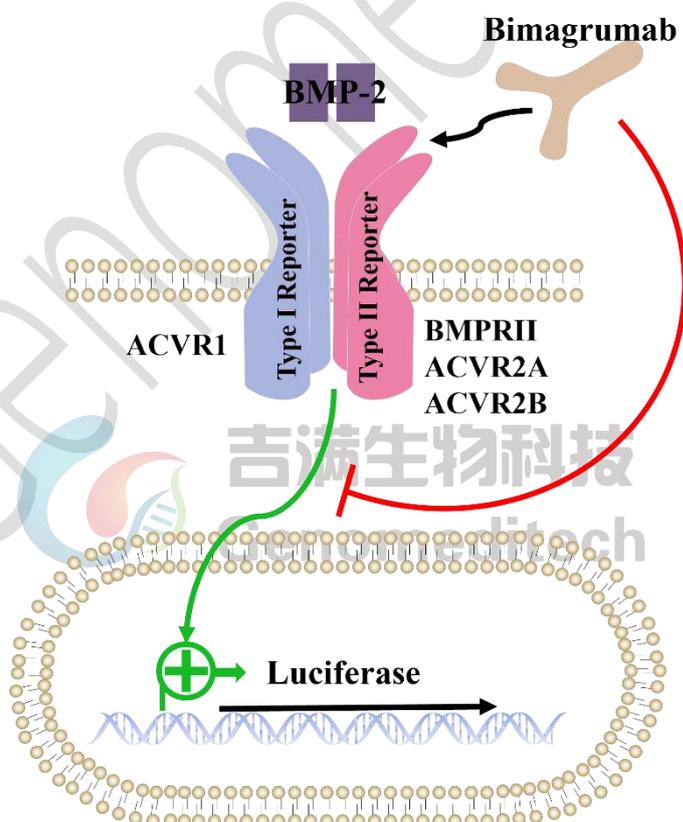

Specifications

Quantity	5E6 Cells per vial, 1 mL
Product Format	1 vial of frozen cells
Shipping	Shipped on dry ice
Storage Conditions	Liquid nitrogen immediately upon receipt

Recovery Medium	DMEM+10% FBS+1% P.S
Growth medium	DMEM+10% FBS+1% P.S+4 µg/mL Blasticidin
Note	None
Freezing Medium	90% FBS+10% DMSO
Growth properties	Adherent
Growth Conditions	37°C, 5% CO ₂

Mycoplasma Testing	The cell line has been screened to confirm the absence of Mycoplasma species.
Safety considerations	Biosafety Level 2
Note	It is recommended to expand the cell culture and store a minimum of 10 vials at an early passage for potential future use.

Materials

Reagent	Manufacturer/Catalogue No.
DMEM	Gibco/C11995500BT
Fetal Bovine Serum	Cegrogen biotech/A0500-3010
Pen/Strep	Thermo/15140-122
Blasticidin	Genomeditech/ GM-040404
Recombinant Human, Mouse, Rat, Rhesus, Canine BMP-2 Protein	Sino Biological/10426-HNAE
Human GDF-2 Protein, Tag Free	ACRO/GD2-H5211
Human Activin A Protein; His Tag	Genomeditech/ GM-87616RP
Human latent GDF-8 Protein; His Tag	Genomeditech/ GM-87733RP
Anti-ACVR2B hIgG1 Antibody(Bimagrumab)	Genomeditech/ GM-51148AB
GMOne-Step Luciferase Reporter Gene Assay Kit	Genomeditech/ GM-040503

Figures

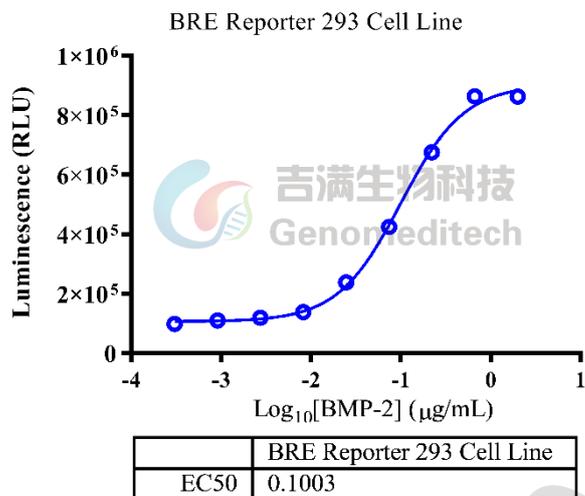


Figure 1 | Response to Recombinant Human, Mouse, Rat, Rhesus, Canine BMP-2 Protein. BRE Reporter 293 Cell Line (Cat. GM-C37665) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Recombinant BMP-2 Protein (Sino Biological/10426-HNAE) in assay buffer (DMEM+1% FBS+1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. [GM-040503](#)). The maximum induction fold was approximately [8.9]. Data are shown by drug mass concentration.

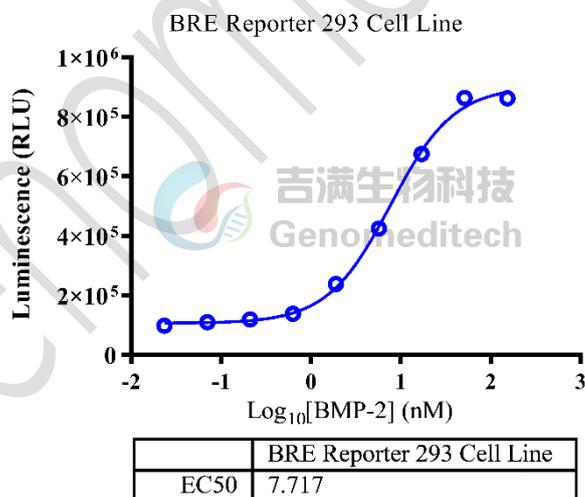


Figure 2 | Response to Recombinant Human, Mouse, Rat, Rhesus, Canine BMP-2 Protein. BRE Reporter 293 Cell Line (Cat. GM-C37665) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Recombinant BMP-2 Protein (Sino Biological/10426-HNAE) in assay buffer (DMEM+1% FBS+1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. [GM-040503](#)). The maximum induction fold was approximately [8.9]. Data are shown by drug molar concentration.

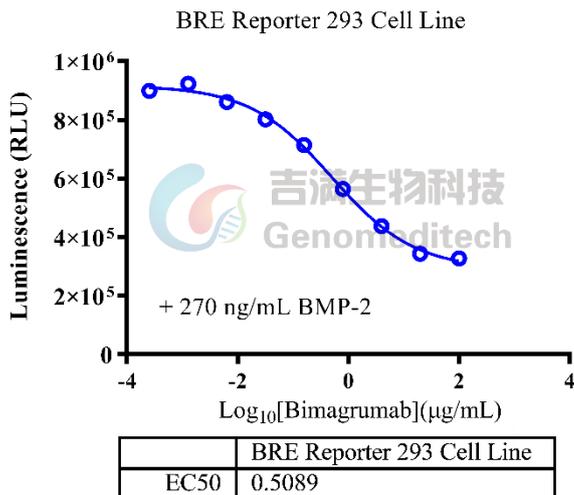


Figure 3 | Response to Anti-ACVR2B hIgG1 Antibody(Bimagrumb). Serial dilutions of the Anti-ACVR2B hIgG1 Antibody(Bimagrumb) (Cat. [GM-51148AB](#)) was incubated with 1.5E4 cells/well of the BRE Reporter 293 Cell Line (Cat. GM-C37665) in a 96-well plate for 1 hour in assay buffer (DMEM + 1% FBS + 1% P.S). Subsequently, the BMP-2 Protein (Sino Biological/10426-HNAE) at a concentration of 27 ng/well was added, and the coculture proceeded for an additional 6 hours. Firefly luciferase activity is then measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. [GM-040503](#)). The results indicated maximum blocking folds of approximately [2.9]. Data are shown by drug mass concentration.

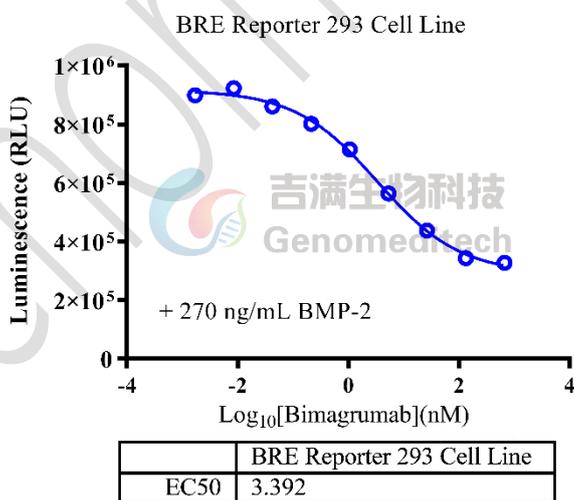


Figure 4 | Response to Anti-ACVR2B hIgG1 Antibody(Bimagrumb). Serial dilutions of the Anti-ACVR2B hIgG1 Antibody(Bimagrumb) (Cat. [GM-51148AB](#)) was incubated with 1.5E4 cells/well of the BRE Reporter 293 Cell Line (Cat. GM-C37665) in a 96-well plate for 1 hour in assay buffer (DMEM + 1% FBS + 1% P.S). Subsequently, the BMP-2 Protein (Sino Biological/10426-HNAE) at a concentration of 27 ng/well was added, and the coculture proceeded for an additional 6 hours. Firefly luciferase activity is then measured using the GMOne-Step Luciferase Reporter

Gene Assay Kit (Cat. [GM-040503](#)). The results indicated maximum blocking folds of approximately [2.9]. Data are shown by drug mmolar concentration.

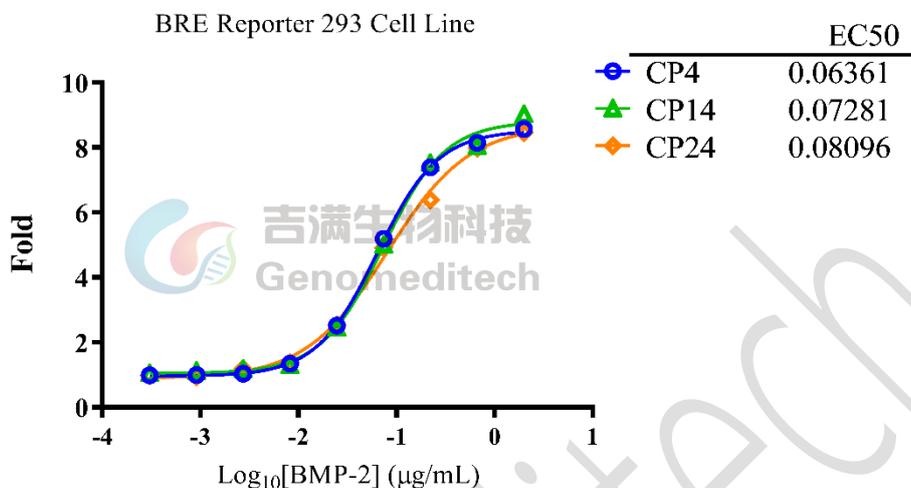


Figure 5 | The passage stability of response to Recombinant BMP-2 Protein. The passage 4, 14, and 24 of BRE Reporter 293 Cell Line (Cat. GM-C37665) at a concentration of 1.5E4 cells/well (96-well format) were stimulated with serial dilutions of Recombinant BMP-2 Protein (Sino Biological/10426-HNAE) in assay buffer (DMEM+1% FBS+1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. [GM-040503](#)). Data are shown by drug mass concentration.

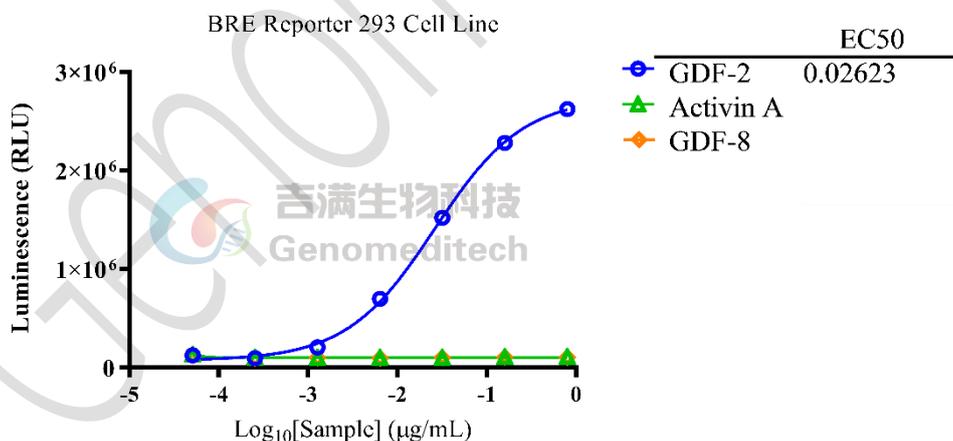


Figure 6 | Response to ligand protein. The BRE Reporter 293 Cell Line (Cat. GGM-C37665) at a concentration of 1.5E4 cells/well (96-well format) was stimulated with serial dilutions of Human GDF-2 Protein (ACRO/GD2-H5211), Human Activin A Protein (Cat. [GM-87616RP](#)) and Human latent GDF-8 Protein(Cat. [GM-87733RP](#)) in assay buffer (DMEM+1% FBS+1% P.S) for 6 hours. The firefly luciferase activity was measured using the GMOne-Step Luciferase Reporter Gene Assay Kit (Cat. [GM-040503](#)). Data are shown by drug mass concentration.

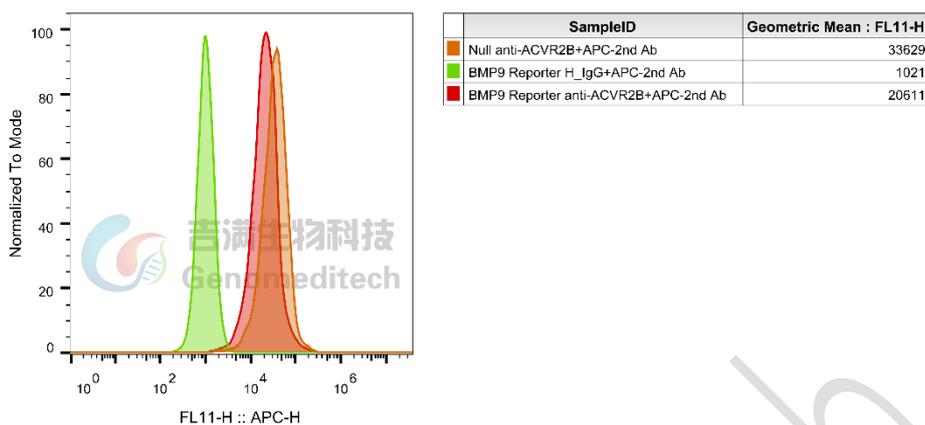


Figure 7 | BRE Reporter 293 Cell Line (Cat. GM-C37665) was determined by flow cytometry using Anti-ACVR2B hIgG1 Antibody(Bimagrumab) (Cat. [GM-51148AB](#)).

Cell Recovery

Recovery Medium: DMEM+10% FBS+1% P.S

To insure the highest level of viability, thaw the vial and initiate the culture as soon as possible upon receipt. If upon arrival, continued storage of the frozen culture is necessary, it should be stored in liquid nitrogen vapor phase and not at -70°C. Storage at -70°C will result in loss of viability.

- Thaw the vial by gentle agitation in a 37°C water bath. To reduce the possibility of contamination, keep the O-ring and cap out of the water. Thawing should be rapid (approximately 2 - 3 minutes).
- Remove the vial from the water bath as soon as the contents are thawed, and decontaminate by dipping in or spraying with 70% ethanol. All of the operations from this point on should be carried out under strict aseptic conditions.
- Transfer the vial contents to a centrifuge tube containing 5.0 mL complete culture medium and spin at approximately 176 x g for 5 minutes. Discard supernatant.
- Resuspend cell pellet with the recommended recovery medium. And dispense into appropriate culture dishes.
- Incubate the culture at 37°C in a suitable incubator. A 5% CO₂ in air atmosphere is recommended if using the medium described on this product sheet.

Cell Freezing

Freezing Medium: 90% FBS+10% DMSO

- Centrifuge at 176 x g for 3 minutes to collect cells.
- Resuspend the cells in pre-cooled freezing medium and adjust the cell density to 5E6 cells/mL.
- Aliquot 1 mL into each vial.
- Place the vial in a controlled-rate freezing container and store at -80°C for at least 1 day, then transfer to liquid nitrogen as soon as possible.

Cell passage

Growth medium: DMEM+10% FBS+1% P.S+4 µg/mL Blasticidin

For the first 1 to 2 passages post-resuscitation, use the recovery medium. Once the cells have stabilized, switch to a growth medium.

- a) Subculturing is necessary when the cell density reaches 80%. It is recommended to perform subculturing at a ratio of 1:3 to 1:4 every 2-3 days. Ensure that the density does not exceed 80%, as overcrowding can lead to reduced viability due to compression.
- b) Remove and discard culture medium.
- c) Briefly rinse the cell layer with PBS to remove all traces of serum that contains trypsin inhibitor.
- d) Add 1.0 mL of 0.25% (w/v) Trypsin-EDTA solution to dish and observe cells under an inverted microscope until cell layer is dispersed (usually within 30 to 60 seconds at 37°C).
- e) Note: To avoid clumping do not agitate the cells by hitting or shaking the flask while waiting for the cells to detach. Cells that are difficult to detach may be placed at 37°C to facilitate dispersal.
- f) Add 2.0 mL of growth medium to mix well and aspirate cells by gently pipetting.
- g) After centrifugation, resuspend the pellet and add appropriate aliquots of the cell suspension to new culture vessels.
- h) Incubate cultures at 37°C.

Subcultivation Ratio: A subcultivation ratio of 1:3 - 1:4 is recommended

Medium Renewal: Every 2 to 3 days

Notes

- a) Upon initial thawing, a higher number of dead cells is observed, which is a normal phenomenon. Significant improvement is seen after adaptation. Once the cells reach a stable state, the number of dead cells decreases after subculturing and the cell growth rate becomes stable.
- b) Ensure that the cell density does not exceed 80%, as overcrowding may lead to reduced viability due to compression.

Related Products

GCGR	
H_GCGR Reporter CHO-K1 Cell Line	H_GCGR Reporter HEK-293 Cell Line
H_GCGR Reporter HEK-293 DDX35TM Cell Line	H_GCGR CHO-K1 Cell Line
H_GCGR HEK-293 Cell Line	Mouse_GCGR HEK-293 Cell Line
Anti-H_GCGR hIgG2 Antibody(volagidemab)	
GLP1R	
H_GLP1R Reporter CHO-K1 Cell Line	H_GLP1R Reporter HEK-293 Cell Line
H_GLP1R Reporter HEK-293 DDX35TM Cell Line	Tango-H_GLP1R CHO-K1 Cell Line
Cynomolgus_GLP1R HEK-293 Cell Line	H_GLP1R CHO-K1 Cell Line
H_GLP1R HEK-293 Cell Line	Mouse_GLP1R HEK-293 Cell Line
Anti-GLP1R hIgG1 Antibody(mAb-36986)	Anti-H_GLP1R hIgG1 Antibody(glutazumab)

FGF21:FGFR	
H_FGF21 Reporter HEK-293 Cell Line	
CALCA(CGRP): CALCRL RAMP	
H_CALCRL RAMP1 Reporter HEK-293 Cell Line	H_CALCRL RAMP1 Reporter HEK-293 DDX35TM Cell Line
Cynomolgus_CALCRL RAMP1 HEK-293 Cell Line	H_CALCRL RAMP1 CHO-K1 Cell Line
H_CALCRL RAMP1 HEK-293 Cell Line	
Anti-CALCRL RAMP1 hIgG2 Antibody(Erenumab)	
GIP:GIPR	
H_GIPR Reporter CHO-K1 Cell Line	H_GIPR Reporter HEK-293 Cell Line
H_GIPR Reporter HEK-293 DDX35TM Cell Line	Cynomolgus_GIPR HEK-293 Cell Line
H_GIPR CHO-K1 Cell Line	H_GIPR HEK-293 Cell Line
Mouse_GIPR HEK-293 Cell Line	
Anti-H_GIPR hIgG1 Antibody(AMG-133)	
ACVR2A: ACTRIIB: Active A	
ACVR2A KO HEK-293 Cell Line	Activin A Reporter Cell Line
H_ACVR2A Reporter Cell Line	H_ACVR2B Reporter Cell Line
ACVR2B KO HEK-293 Cell Line	H_ACVR2A HEK-293(ACVR2B KO) Cell Line
H_ACVR2B CHO-K1 Cell Line	H_ACVR2B HEK-293(ACVR2A KO) Cell Line
Anti-ACVR2B hIgG1 Antibody(Bimagrumab)	Anti-ACVR2B hIgG1 Antibody(Fab-17G05)
Anti-ACVR2B mIgG2a Antibody(Bimagrumab)	Anti-H_ACVR2B hIgG1 Reference Antibody(Bimbio)
Biotinylated Human ACVR2A Protein; His-Avi Tag	Biotinylated Human ACVR2B Protein; His-Avi Tag
Biotinylated Mouse ACVR2A Protein; His-Avi Tag	Biotinylated Mouse ACVR2B Protein; His-Avi Tag
Human Activin A Protein; His Tag	Human Activin B Protein; His Tag
Human ACVR2A Protein; hFc Tag	Human ACVR2A Protein; His Tag
Human ACVR2B Protein; hFc Tag	Human ACVR2B Protein; His Tag
Human latent GDF-8 Protein; His Tag	Mouse ACVR2B Protein; His Tag
AMY: CALCR RAMP	
H_CALCR RAMP3(AMY3) Reporter CHO-K1 Cell Line	H_CALCR Reporter CHO-K1 Cell Line

Limited Use License Agreement

Genomeditech (Shanghai) Co., Ltd grants to the Licensee all intellectual property rights, exclusive, non-transferable, and non-sublicensable rights of the Licensed Materials; Genomeditech (Shanghai) Co., Ltd will retain ownership of the Licensed Materials, cell line history packages, progeny, and the Licensed Materials including modified materials.

Between Genomeditech (Shanghai) Co., Ltd, and Licensee, Licensee is not permitted to modify cell lines in any way. The Licensee shall not share, distribute, sell, sublicense, or otherwise provide the Licensed Materials, or progenitors to third parties such as laboratories, departments, research institutions, hospitals, universities, or biotechnology companies for use other than for the purpose of outsourcing the Licensee's research.

Please refer to the Genomeditech Cell Line License Agreement for details.

Genomeditech